



BD Infinicyt™ Software

Advanced flow cytometry analysis software
designed for the way you work



BD Infinicyt™ Software

Advanced analysis software with innovative tools for flow cytometry data analysis



Complex Analysis



Standardization and Efficiency



Safety and Security




BD Total Research Solution

Take your research to the next level with advanced features

BD Infinicyt™ Software offers a variety of classic and multiparametric diagrams, analysis profiles [templates], analysis strategies and statistical outputs. Additional key advanced features include:


- Automatic Population Separator (APS)

Based on Principal Component Analysis (PCA), APS uses the information from all the parameters in the file to expose the cell populations with the best separation.

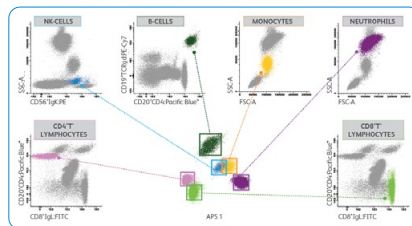


- Clustering

Patented algorithms that automatically identify similar cell populations and group them into clusters for further exploration.

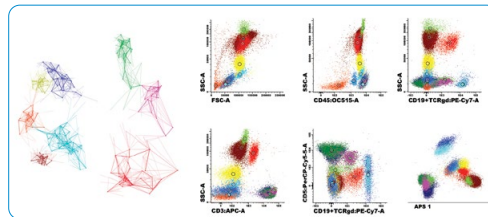

- Maturation patterns

Define maturation pathways in your research samples (e.g. myeloid development research). Allow the software to automatically adjust and divide the events into population stages.

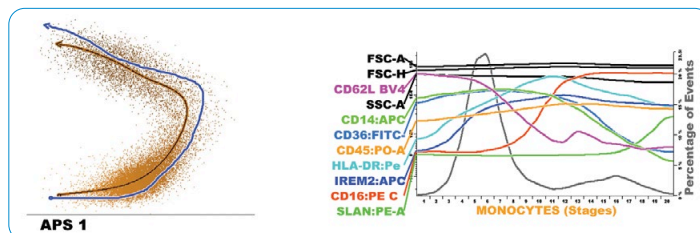

- File merging

Combine multiple files, from the same or different samples, and visualize them all on the workspace. All files can be analyzed at the same time.

- Reference images

Visualize and compare your samples to saved references (e.g., a negative control or a normal population), or track treated samples over time.

Next Generation Flow™ Solution

Automate your analysis and increase standardization and efficiency¹ when using predefined assays

The Next Generation Flow™ Solution is based on the standardization of your entire workflow, from sample preparation of predefined panels to expert-guided automatic reports. BD Infinicyt™ Software provides automated analysis and reporting support for Next Generation Flow™ research panels.



- Off-the-shelf research assay panels

Templates and databases for: PIDOT (Primary Immunodeficiency Orientation Tube), ALOT (Acute Leukemia Orientation Tube), LST (Lymphoid Screening Tube), BCP-ALL MRD (B-Cell Precursor Acute Lymphocytic Leukemia Minimal Residual Disease), MM MRD (Multiple Myeloma Minimal Residual Disease), and CPC (Circulating Plasma Cells).

- Automated gating and identification

Helps to identify normal cell populations by comparing to databases built with normal reference samples. Automated clustering of events and immune profiling of samples.

- Automatic Reports

Generates comments and conclusions based on results compared to internal database reference ranges. Immunophenotype of the abnormal population provided. Sample-specific Limit of Detection (LOD) and Lower Limit of Quantitation (LLOQ) definitions and sample quality assessments (e.g., hemodilution) provided.

Tools for your compliance and for product safety and security:

- Features supporting 21 CFR part 11 compliance
- User-access control and authentication of unique users
- Password-protected analyzed files
- Audit trails
- Electronic signatures
- Anonymization and encryption
- Auto-recover tool

Add to the BD total research solution with BD Infinicyt™ Software:

Simplify your procurement processes and receive world-class service and support on your entire workflow.

- End-to-end experiment capabilities from a single vendor
Reagents + Sample Preparation + Flow Cytometers + Advanced Analysis Software
- Worldwide support
Capabilities to troubleshoot through many components of your flow cytometry workflow.
- 45+ years of experience in flow cytometry analysis
Support your advanced application needs

BD Infinicyt™ Software licensing options are available to meet your configuration needs. Reach out to a BD representative.

Licensing options are available

- User formats: [Standalone and network](#)
- Security formats: [Hardware \(dongle\) and software activation key](#)



Visit our website for more information
bdbiosciences.com/infinicyt

References

- 1 Flores-Montero J, Sanoja-Flores L, Paiva B, et al. Next Generation Flow for highly sensitive and standardized detection of minimal residual disease in multiple myeloma. *Leukemia*. 2017;31(10):2094-2103. doi: 10.1038/leu.2017.29

For Research Use Only. Not for use in diagnostic or therapeutic procedures.

BD Life Sciences, Milpitas, CA 95035, U.S.

bdbiosciences.com

BD, the BD Logo, Infinicyt and Next Generation Flow are trademarks of Becton, Dickinson and Company or its affiliates. All other trademarks are the property of their respective owners. © 2023 BD. All rights reserved. BD-93620 (v1.0) 0923

