

Technical Data Sheet


FITC Rabbit Anti- Active Caspase-3

Product Information

Material Number:	559341
Alternate Name:	CPP32; Yama; Apopain
Size:	100 tests
Vol. per Test:	20 µl
Clone:	C92-605
Isotype:	Rabbit IgG
Reactivity:	QC Testing: Human Tested in Development: Mouse
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The caspase family of cysteine proteases plays a key role in apoptosis and inflammation. Caspase-3 is a key protease that is activated during the early stages of apoptosis and, like other members of the caspase family, is synthesized as an inactive pro-enzyme that is processed in cells undergoing apoptosis by self-proteolysis and/or cleavage by another protease. The processed forms of caspases consist of large (17-22 kDa) and small (10-12 kDa) subunits which associate to form an active enzyme. Active caspase-3, a marker for cells undergoing apoptosis, consists of a heterodimer of 17 and 12 kDa subunits which is derived from the 32 kDa pro-enzyme. Active caspase-3 proteolytically cleaves and activates other caspases, as well as relevant targets in the cytoplasm, e.g., D4-GDI and Bcl-2, and in the nucleus (e.g. PARP). This antibody has been reported to specifically recognize the active form of caspase-3 in human and mouse cells. It has not been reported to recognize the pro-enzyme form of caspase-3.

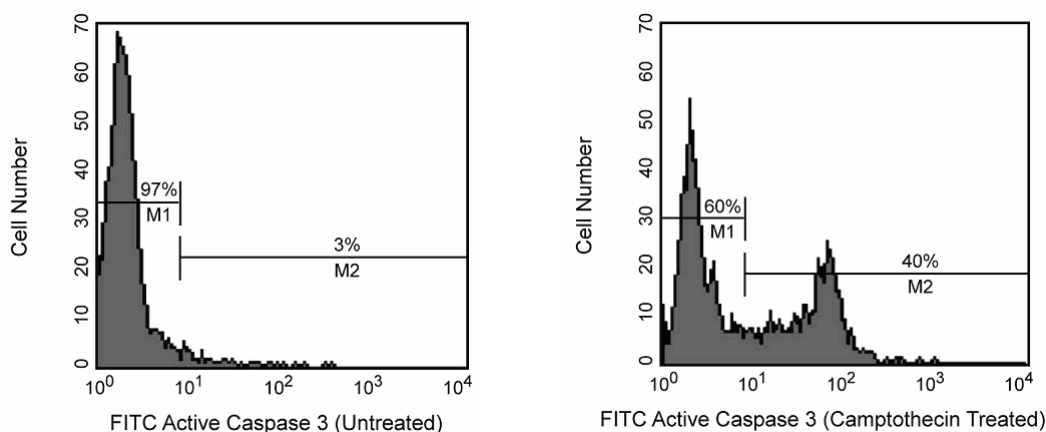

Flow cytometric analysis of apoptotic and non-apoptotic populations for active caspase-3. Jurkat cells (Human T-cell leukemia; ATCC TIB-152) were left untreated (left panel) or treated with 4 µM of camptothecin for 4 hr to induce apoptosis (right panel). Cells were washed once in PBS, then fixed and permeabilized using the BD Cytofix/Cytoperm™ Kit (Cat. No. 554714) for 20 min at room temperature (RT), pelleted and washed with BD Perm/Wash™ buffer (component of Cat. No. 554714). Cells were subsequently stained with the FITC rabbit anti- active caspase-3 antibody (clone C92-605). Cells were then washed and resuspended in BD Perm/Wash™ buffer before analyzing by flow cytometry. The results show that untreated cells were primarily negative for active caspase-3 (left panel, M1); whereas over one third of the treated cells were positive for active caspase-3 staining (right panel, M2).

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



Application Notes

Application

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

Suggested Companion Products

<u>Catalog Number</u>	<u>Name</u>	<u>Size</u>	<u>Clone</u>
554714	BD Cytotfix/Cytoperm™ Fixation/Permeablization Kit	250 tests	(none)

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
3. For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/pharmingen/colors.
4. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

References

- Dai C, Krantz SB. Interferon gamma induces upregulation and activation of caspases 1, 3, and 8 to produce apoptosis in human erythroid progenitor cells. *Blood*. 1999; 93(10):3309-3316.(Biology)
- Fujita N, Tsuruo T. Involvement of Bcl-2 cleavage in the acceleration of VP-16-induced U937 cell apoptosis. *Biochem Biophys Res Commun*. 1998; 246(2):484-488.(Biology)
- Thornberry NA, Lazebnik Y. Caspases: enemies within. *Science*. 1998; 281(5381):1312-1316.(Biology)