

Technical Data Sheet


FITC Mouse Anti-Human β 2-Microglobulin**Product Information**

Material Number:	551338
Size:	100 Tests
Vol. per Test:	20 μ l
Clone:	TÜ99
Isotype:	Mouse IgM, κ
Reactivity:	QC Testing: Human
Workshop:	N/A
Storage Buffer:	Aqueous buffered solution containing BSA and $\leq 0.09\%$ sodium azide.

Description

β 2-microglobulin is the 12 kDa β chain which noncovalently associates with the 44 kDa α chain to form the HLA Class I antigen complex. This antibody reacts with free or complexed β 2-microglobulin. TÜ99 does not react with β 2-microglobulin of other species. This antibody is suitable for staining formalin-fixed, paraffin-embedded tissue sections using TUF pretreatment.

This antibody is routinely tested by flow cytometric analysis. Other applications were tested at BD Biosciences Pharmingen during antibody development only or reported in the literature.

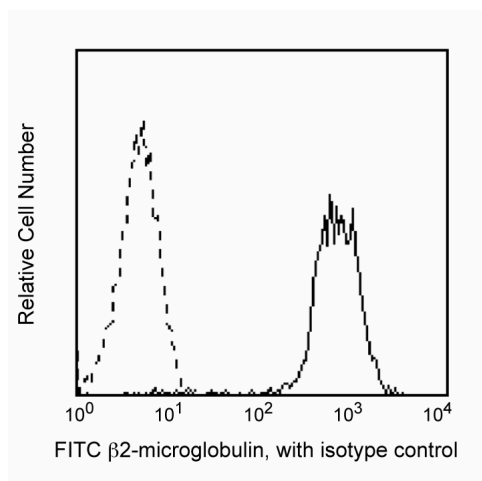

Profile of peripheral blood lymphocytes analyzed by flow cytometry

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed. Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes**Application**

Flow cytometry	Routinely Tested
----------------	------------------

Suggested Companion Products

<u>Catalog Number</u>	<u>Name</u>	<u>Size</u>	<u>Clone</u>
555583	FITC Mouse IgM, κ Isotype Control	100 Tests	G155-228

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
3. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
© 2017 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.

551338 Rev. 4



4. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.

References

- Barclay NA, Brown MH, Birkeland ML, et al, ed. *The Leukocyte Antigen FactsBook*. San Diego, CA: Academic Press; 1997(Biology)
- Desoye G, Dohr GA, Motter W, et al. Lack of HLA class I and class II antigens on human preimplantation embryos. *J Immunol*. 1988; 140(12):4157-4159. (Biology)
- Dohr GA, Motter W, Leitinger S, et al. Lack of expression of HLA [corrected] class I and class II molecules on the human oocyte. *J Immunol*. 1987; 138(11):3766-3770. (Biology)