

Technical Data Sheet


APC Mouse Anti-Human CD25

Product Information

Material Number:	555434
Alternate Name:	IL-2R; IL2RA; IL-2R α ; TCGFR; TAC antigen; p55
Size:	100 Tests
Vol. per Test:	20 μ l
Clone:	M-A251
Immunogen:	Phytohemagglutinin stimulated human lymphocytes
Isotype:	Mouse (BALB/c) IgG1, κ
Reactivity:	QC Testing: Human Tested in Development: Rhesus, Cynomolgus, Baboon
Workshop:	IV A053
Storage Buffer:	Aqueous buffered solution containing BSA and \leq 0.09% sodium azide.

Description

The M-A251 monoclonal antibody specifically binds to the 55 kDa type I transmembrane glycoprotein known as low-affinity interleukin-2 receptor alpha chain subunit (IL-2R α). CD25 is expressed on regulatory T cells, activated lymphocytes (T and B), and monocytes. It associates with the IL-2R β /CD122 and IL-2R γ /CD132 receptor chains to form the high-affinity IL-2R complex. CD25 expression on T and B lymphocytes is upregulated by antigenic or mitogenic stimulation. Soluble CD25/IL-2R α is produced as a consequence of lymphocyte stimulation and is found in biological fluids following inflammatory responses.

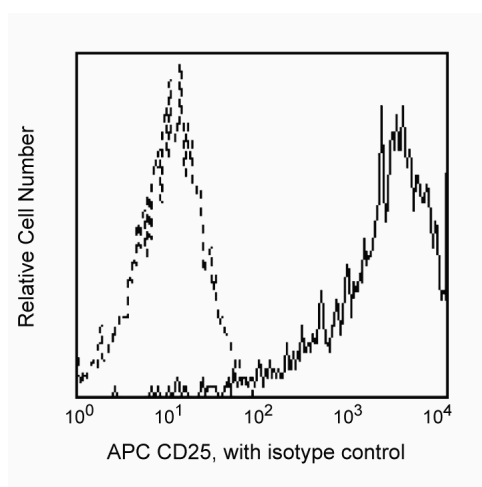

Flow cytometric analysis of CD25 expression on PHA-stimulated human peripheral blood lymphocytes.
Phytohemagglutinin-stimulated (3 days) PBMCs were stained with either APC Mouse Anti-Human CD25 (Cat. No. 555434/561399/560987; solid line histogram) or APC Mouse IgG1 κ Isotype Control (Cat. No. 555751; dashed line histogram). Fluorescent histograms were derived from gated events with the side and forward light-scattering characteristics of viable lymphocytes. Flow cytometry was carried out on a BD FACScan™ system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

Application Notes

Application

Flow cytometry	Routinely Tested
----------------	------------------

Suggested Companion Products

Catalog Number	Name	Size	Clone
555751	APC Mouse IgG1, κ Isotype Control	100 Tests	MOPC-21
561399	APC Mouse Anti-Human CD25	50 Tests	M-A251
554656	Stain Buffer (FBS)	500 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)
560987	APC Mouse Anti-Human CD25	25 Tests	M-A251

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
© 2017 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.

555434 Rev. 9



Product Notices

1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
5. This APC-conjugated reagent can be used in any flow cytometer equipped with a dye, HeNe, or red diode laser.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
7. Species testing during development may have been performed with a different format of the same clone. Selected applications have been tested for cross-reactivity.
8. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

Knapp W. W. Knapp .. et al., ed. *Leucocyte typing IV : white cell differentiation antigens*. Oxford New York: Oxford University Press; 1989:1-1182. (Clone-specific)
Schlossman SF. Stuart F. Schlossman .. et al., ed. *Leucocyte typing V : white cell differentiation antigens : proceedings of the fifth international workshop and conference held in Boston, USA, 3-7 November, 1993*. Oxford: Oxford University Press; 1995(Biology)