

## Technical Data Sheet


## APC Mouse Anti-Human CD184

## Product Information

<b>Material Number:</b>	555976
<b>Alternate Name:</b>	CXCR4; Fusin; SDF-1 receptor; LAP3; LCR1; LESTR; NPY3R; NPY3R; WHIM; HM8
<b>Size:</b>	100 tests
<b>Vol. per Test:</b>	20 µl
<b>Clone:</b>	12G5
<b>Isotype:</b>	Mouse (BALB/c) IgG2a, κ
<b>Reactivity:</b>	QC Testing: Human
<b>Workshop:</b>	VII 70204, 70305
<b>Storage Buffer:</b>	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

## Description

The 12G5 monoclonal antibody specifically binds to CD184, also known as CXCR4 and Fusin. CD184/CXCR4 is a seven-transmembrane domain, G-protein-linked, glycoprotein chemokine receptor. CD184 serves as a receptor for the C-X-C chemokine, SDF-1. It is expressed on a wide variety of hematopoietic cells, vascular endothelial cells and cells of the nervous system. CD184 plays a variety of roles in hematopoiesis, vascularization and neural development. CD184 also functions as a coreceptor for infection with T-cell tropic strains of HIV-1 and as a receptor for CD4-independent infection by some HIV isolates. The 12G5 antibody has been reported to block CD4-independent infection by HIV-2 and CD4-dependent infection by some T-cell tropic isolates of HIV-1.

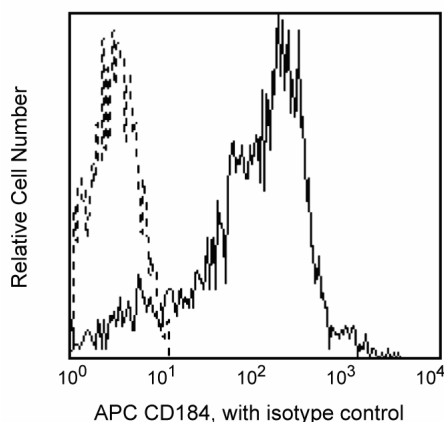

Profile of peripheral blood lymphocytes analyzed by flow cytometry.

## Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to APC under optimum conditions, and unconjugated antibody and free APC were removed.

## Application Notes

## Application

Flow cytometry	Routinely Tested
----------------	------------------

## Recommended Assay Procedure:

**Flow cytometry:** Immunophenotyping studies of chemokine receptors need to be performed on freshly collected whole blood (<24 hours). Incubation with the antibody should be done at room temperature in the dark. Cellular manipulation, such as Ficoll separation, freezing, or exposure to cold temperatures prior to staining have been shown to cause a decrease in staining intensity and inconsistent results.

## Suggested Companion Products

Catalog Number	Name	Size	Clone
555576	APC Mouse IgG2a, κ Isotype Control	100 tests	G155-178

## BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



## Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100- $\mu$ l experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
5. This APC-conjugated reagent can be used in any flow cytometer equipped with a dye, HeNe, or red diode laser.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
7. Please refer to [www.bdbiosciences.com/pharming/protocols](http://www.bdbiosciences.com/pharming/protocols) for technical protocols.

## References

- Endres MJ, Clapham PR, Marsh M, et al. CD4-independent infection by HIV-2 is mediated by fusin/CXCR4. *Cell*. 1996; 87(4):745-756. (Biology)
- Feng Y, Broder CC, Kennedy PE, Berger EA. HIV-1 entry cofactor: functional cDNA cloning of a seven-transmembrane, G protein-coupled receptor. *Science*. 1996; 272(5263):872-877. (Biology)
- Loetscher M, Geiser T, O'Reilly T, Zwahlen R, Baggiolini M, Moser B. Cloning of a human seven-transmembrane domain receptor, LESTR, that is highly expressed in leukocytes. *J Biol Chem*. 1994; 269(1):232-237. (Biology)
- Nagasawa T, Nakajima T, Tachibana K, et al. Molecular cloning and characterization of a murine pre-B-cell growth-stimulating factor/stromal cell-derived factor 1 receptor, a murine homolog of the human immunodeficiency virus 1 entry coreceptor fusin. *Proc Natl Acad Sci U S A*. 1996; 93(25):14726-14729. (Biology)
- Simmons G, Wilkinson D, Reeves JD, et al. Primary, syncytium-inducing human immunodeficiency virus type 1 isolates are dual-tropic and most can use either Lestr or CCR5 as coreceptors for virus entry. *J Virol*. 1996; 70(12):8355-8360. (Biology)

## BD Biosciences

[bdbiosciences.com](http://bdbiosciences.com)

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.268.5430	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD

