

Technical Data Sheet


Alexa Fluor® 700 Mouse Anti-Human CD4

Product Information

Material Number:	557922
Size:	0.1 mg
Concentration:	0.2 mg/ml
Clone:	RPA-T4
Isotype:	Mouse IgG1, κ
Reactivity:	QC Testing: Human
Workshop:	IV T114
Storage Buffer:	Aqueous buffered solution containing protein stabilizer and ≤0.09% sodium azide.

Description

The RPA-T4 clone reacts with CD4, a 59 kDa single-chain transmembrane glycoprotein [receptor for human immunodeficiency virus (HIV)] present on T-helper/inducer cell populations. This antibody binds to the D1 domain (CDR1 and CDR3 epitopes) of the CD4 antigen and reacts with approximately 80% of thymocytes and 45% of peripheral blood lymphocytes. CD4 is also present in low density on peripheral blood monocytes. RPA-T4 is capable of blocking HIV-1, gp120, and inhibits syncytium formation. This antibody is routinely tested by flow cytometric analysis. Other applications were tested at BD Biosciences Pharmingen during antibody development only or reported in the literature.

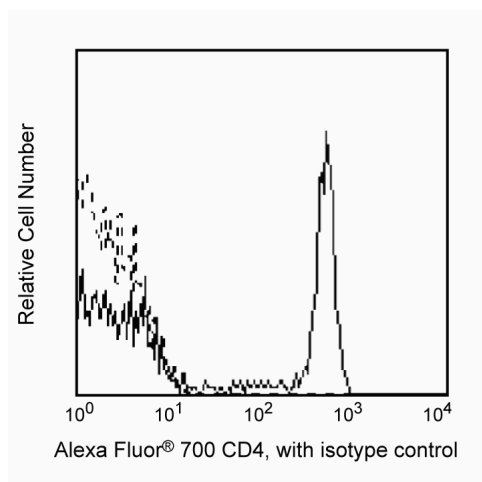

Profile of CD4 (RPA-T4) reactivity on peripheral blood lymphocytes analyzed by flow cytometry.

Preparation and Storage

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 700 under optimum conditions, and unreacted Alexa Fluor® 700 was removed.

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

Application Notes

Application

Flow cytometry	Routinely Tested
----------------	------------------

Suggested Companion Products

Catalog Number	Name	Size	Clone
557882	Alexa Fluor® 700 Mouse IgG1, κ Isotype Control	0.1 mg	MOPC-21

Product Notices

- This antibody has been developed for the application listed above. However, a routine test is not performed on every lot. Researchers should determine the optimal concentration of this reagent for their individual applications.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- For fluorochrome spectra and suitable instrument settings, please refer to our Fluorochrome Web Page at www.bdbiosciences.com/colors.

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



4. The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
5. Alexa Fluor® 700 has an adsorption maximum of ~700nm and a peak fluorescence emission of ~720nm. Before staining cells with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
6. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
7. Alexa Fluor is a registered trademark of Molecular Probes, Inc., Eugene, OR.

References

- Beavis AJ, Pennline KJ. Allo-7: a new fluorescent tandem dye for use in flow cytometry. *Cytometry*. 1996; 24(4):390-395. (Biology)
- Knapp W, Dorken B, Rieber EP, et al, ed. *Leukocyte Typing IV*. New York: Oxford University Press; 1989. (Clone-specific)
- Roederer M, Kantor AB, Parks DR, Herzenberg LA. Cy7PE and Cy7APC: bright new probes for immunofluorescence. *Cytometry*. 1996; 24(3):191-197. (Biology)
- Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. New York: Oxford University Press; 1995. (Clone-specific)