

Technical Data Sheet

Alexa Fluor® 647 Rat anti-Histone H3 (pS28)


Product Information

Material Number:	558609
Size:	100 tests
Vol. per Test:	5 µl
Clone:	HTA28
Immunogen:	Phosphorylated Human Histone H3 Peptide
Isotype:	Rat IgG2a, κ
Reactivity:	Tested: Human Reported: Cow, Drosophila, Mouse, Rat
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

Histones are highly basic proteins that complex with DNA to form chromatin. Histone H3 is a ~15-kDa protein that is phosphorylated at serine 28 (S28), S10, and/or threonine 11 during mammalian cell mitosis and meiosis. The phosphorylation sites are located in the N-terminal tail, a region that is outside of the chromatin fiber and is thus accessible for interactions with agents that may regulate chromatin or specific gene activities. The phosphorylation states of the two serine sites during the cell cycle are highly regulated by Aurora B kinase and a PP1 phosphatase: S10 is in the phosphorylated state from late G2 phase to anaphase, while S28 is phosphorylated from prophase to anaphase. Furthermore, phosphorylation of histone H3 S28 may be mediated by other kinases in response to external stimuli. Evidence suggests that histone phosphorylation is involved in the regulation of chromosome condensation, cell division, and gene transcription.

The HTA28 monoclonal antibody reacts with histone H3 phosphorylated at S28 in its N-terminal tail. It does not recognize the unphosphorylated protein.

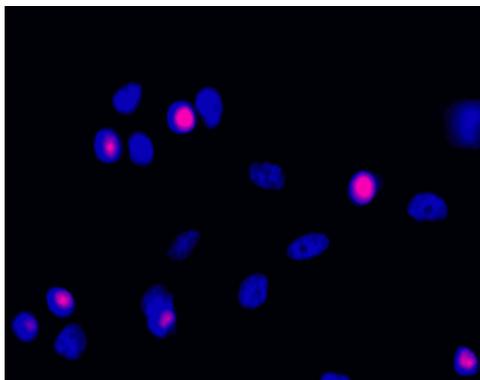

Immunofluorescent staining of human cell lines. HeLa cells (ATCC CCL-2) were cultured overnight and treated for 4 hours with 15 ng/ml colchicine to increase the proportion of mitotic cells. Then they were fixed, permeabilized with cold methanol, stained with Alexa Fluor® 647 Rat anti-Histone H3 (pS28) (pseudo-colored red, which appears pink when co-localized with the blue) and counter-stained with Hoechst 33342 (pseudo-colored blue) according to the Recommended Assay Procedure. The images were captured on a BD Pathway™ 855 Bioimager System with a 20x objective and merged using BD Attovision™ software. This antibody also stains A549 (ATCC CCL-185) and U-2 OS (ATCC HTB-96) cells, and it works with either cold methanol or Triton X-100 permeabilization (see Recommended Assay Procedure).

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated to Alexa Fluor® 647 under optimum conditions, and unreacted Alexa Fluor® 647 was removed.

Application Notes

Application

Bioimaging	Routinely Tested
------------	------------------

Recommended Assay Procedure:

- Seed the cells in appropriate culture medium at ~10,000 cells per well in a BD Falcon™ 96-well Imaging Plate (Cat. No. 353219), and culture overnight.
- Remove the culture medium from the wells, and fix the cells by adding 100 µl of fresh 3.7% Formaldehyde in PBS or BD Cytotfix™ fixation buffer (Cat. No. 554655) to each well and incubating for 10 minutes at room temperature (RT).
- Remove the fixative from the wells, and permeabilize the cells using either cold methanol or Triton™ X-100:
 - Add 100 µl of -20°C 90% methanol or -20°C BD™ Phosflow Perm Buffer III (Cat. No. 558050) to each well and incubate for 5 minutes at RT.

OR

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	888.259.0187	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country-specific contact information, visit bdbiosciences.com/how_to_order/

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

BD, BD Logo and all other trademarks are the property of Becton, Dickinson and Company. ©2008 BD



- b. Add 100 μ l of 0.1% Triton™ X-100 to each well and incubate for 5 minutes at RT.
Triton is a trademark of The Dow Chemical Company.
- Remove the permeabilizer, and wash the wells twice with 100 μ l of 1 \times PBS.
 - Remove the PBS, and block the cells by adding 100 μ l of blocking buffer (3% FBS in 1 \times PBS) or BD Pharmingen™ Stain Buffer (FBS) (Cat. No. 554656) to each well and incubating for 30 minutes at RT.
 - Remove the blocking buffer, dilute the antibody conjugate 1:10 in blocking buffer or Stain Buffer (FBS), and stain the cells by adding 50 μ l of the diluted antibody conjugate to each well and incubating for 1 hour at RT.
 - Remove the diluted antibody conjugate, and wash the wells three times with 100 μ l of 1 \times PBS.
 - Remove the PBS, and counter-stain the nuclei by adding 100 μ l of a 2 μ g/ml solution of Hoechst 33342 (eg, Sigma-Aldrich Cat. No. B2261) in 1 \times PBS to each well at least 15 minutes before imaging.
 - View and analyze the cells on an appropriate imaging instrument.

Suggested Companion Products

Catalog Number	Name	Size	Clone
353219	BD Falcon™ 96-well Imaging Plate	1 box	(none)
554656	Stain Buffer (FBS)	500 ml	(none)
554655	Fixation Buffer	100 ml	(none)
558050	Perm Buffer III	125 ml	(none)

Product Notices

- This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
- This product is sold under license from Shigei Medical Research Institute, Okayama, Japan.
- Alexa Fluor is a registered trademark of Molecular Probes, Inc., Eugene, OR.
- The Alexa Fluor®, Pacific Blue™, and Cascade Blue® dye antibody conjugates in this product are sold under license from Molecular Probes, Inc. for research use only, excluding use in combination with microarrays, or as analyte specific reagents. The Alexa Fluor® dyes (except for Alexa Fluor® 430), Pacific Blue™ dye, and Cascade Blue® dye are covered by pending and issued patents.
- Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- Please refer to wwwbdbiosciences.com/pharmingen/protocols for technical protocols.

References

- Carmena M, Earnshaw WC. The cellular geography of aurora kinases. *Nat Rev Mol Cell Biol.* 2003; 4:842-854. (Biology)
- Cheung P, Allis CD, Sassone-Corsi P. Signaling to chromatin through histone modifications. *Cell.* 2000; 103:263-271. (Biology)
- Choi HS, Choi BY, Cho Y-Y, Zhu F, Bode AM, Dong Z. Phosphorylation of Ser28 in histone H3 mediated by mixed lineage kinase-like mitogen-activated protein triple kinase α . *J Biol Chem.* 2005; 280(14):13545-13553. (Biology)
- Dyson MH, Thomson S, Inagaki M, et al. MAP kinase-mediated phosphorylation of distinct pools of histone H3 at S10 or S28 via mitogen- and stress-activated kinase 1/2. *J Cell Sci.* 2005; 118(10):2247-2259. (Clone-specific: Immunofluorescence)
- Furukawa K, Sugiyama S, Osouda S, et al. Barrier-to-autointegration factor plays crucial roles in cell cycle progression and nuclear organization in *Drosophila*. *J Cell Sci.* 2003; 116(18):3811-3823. (Clone-specific)
- Goto H, Tomono Y, Ajiro K, et al. Identification of a novel phosphorylation site on histone H3 coupled with mitotic chromosome condensation. *J Biol Chem.* 1999; 274(36):25543-25549. (Immunogen: Immunofluorescence)
- Hirata A, Inada K, Tsukamoto T, et al. Characterization of a monoclonal antibody, HTA28, recognizing a histone H3 phosphorylation site as a useful marker of M-phase cells. *J Histochem Cytochem.* 2004; 52(11):1503-1509. (Clone-specific)
- Nowak SJ, Corces VG. Phosphorylation of histone H3: a balancing act between chromosome condensation and transcriptional activation. *Trends Genet.* 2004; 20(4):214-220. (Biology)
- Pascreau G, Arlot-Bonnemains Y, Prigent C. Phosphorylation of histone and histone-like proteins by aurora kinases during mitosis. *Prog Cell Cycle Res.* 2003; 5:369-374. (Biology)