

Technical Data Sheet

V500 Mouse anti-Human HLA-DR**Product Information**


Material Number:	561225
Alternate Name:	MHC class II antigen; HLA class II histocompatibility antigen
Size:	25 Tests
Vol. per Test:	5 µl
Clone:	G46-6
Isotype:	Mouse IgG2a, κ
Reactivity:	QC Testing: Human Tested in Development: Rhesus, Cynomolgus, Baboon Reported Reactivity: Dog
Storage Buffer:	Aqueous buffered solution containing protein stabilizer, glycerol and ≤0.09% sodium azide.

Description

The G46-6 monoclonal antibody specifically binds to HLA-DR, a major histocompatibility complex (MHC) class II antigen. HLA-DR antigens are encoded by genes within the Human Leukocyte Antigen (HLA) Complex located on chromosome 6. HLA-DR is a transmembrane heterodimeric glycoprotein composed of an α chain (36 kDa) and a β subunit (27 kDa) expressed primarily on antigen presenting cells: B cells, dendritic cells, monocytes, macrophages, and thymic epithelial cells. HLA-DR is also expressed on activated T cells. This molecule plays a major role in mediating cellular interactions during antigen presentation to CD4-positive T cells.

The antibody is conjugated to BD Horizon™ V500, which has been developed for use in multicolor flow cytometry experiments and is available exclusively from BD Biosciences. It is excited by the Violet laser with an Ex max of 415 nm and Em Max at 500 nm. BD Horizon V500 conjugates emit at a similar wavelength to Amcyan yet exhibit reduced spillover into the FITC channel. For more information on BD Horizon V500, visit bdbiosciences.com/colors.

When compensating dyes in this spectral range (such as Horizon™ V500 and AmCyan), the most accurate compensation can be obtained using single stained cellular controls. Due to spectral differences between cells and beads in this channel, using BD CompBeads can result in spillover errors for V500 and AmCyan reagents. Therefore, the use of BD CompBeads or BD CompBeads Plus to determine spillover values for these reagents is not recommended. Different V500 reagents (e.g. CD4 vs. CD45) can have slightly different fluorescence spillover therefore, it may also be necessary to use clone specific compensation controls when using these reagents.

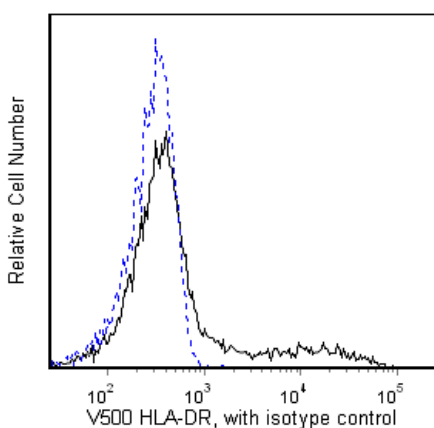

Flow cytometric analysis of HLA-DR expressed on human lymphocytes. Whole blood was stained with BD Horizon™ V500 Mouse anti-Human HLA-DR antibody (Cat. No. 561225/561224; solid line histogram) or with a BD Horizon™ V500 Mouse IgG2a, κ Isotype control (Cat. No. 561221; dotted line histogram). The erythrocytes were lysed with BD PharmLyse™ Lysing Buffer (Cat. No. 555899). The fluorescence histograms were derived from events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometry was performed using a BD FACSCanto™ II Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ V500 under optimum conditions, and unreacted BD Horizon™ V500 was removed.

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
© 2017 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.



Application Notes

Application

Flow cytometry

Routinely Tested

Suggested Companion Products

<u>Catalog Number</u>	<u>Name</u>	<u>Size</u>	<u>Clone</u>
561221	V500 Mouse IgG2a, κ Isotype control	0.1 mg	G155-178
555899	Lysing Buffer	100 mL	(none)
554656	Stain Buffer (FBS)	500 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)
349202	BD FACSTM Lysing Solution	100 mL	(none)
561224	V500 Mouse anti-Human HLA-DR	100 Tests	G46-6

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100- μ l experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. BD Horizon V500 has a maximum absorption of 415 nm and maximum emission of 500 nm. Before staining with this reagent, please confirm that your flow cytometer is capable of exciting the fluorochrome and discriminating the resulting fluorescence.
5. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
6. Species testing during development may have been performed with a different format of the same clone. Selected applications have been tested for cross-reactivity.
7. BD Horizon V500 is covered by one or more of the following US patents: 8,431,416.
8. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.

References

Barclay NA, Brown MH, Birkeland ML, et al, ed. *The Leukocyte Antigen FactsBook*. San Diego, CA: Academic Press; 1997(Biology)

Dieckmann D, Plottner H, Berchtold S, Berger T, Schuler G. Ex vivo isolation and characterization of CD4(+)CD25(+) T cells with regulatory properties from human blood. *J Exp Med*. 2001; 193(11):1303-1310. (Clone-specific: Flow cytometry)

Ibisch C, Pradal G, Bach JM, Lieubeau B. Functional canine dendritic cells can be generated in vitro from peripheral blood mononuclear cells and contain a cytoplasmic ultrastructural marker. *J Immunol Methods*. 2005; 298(1-2):175-82. (Clone-specific)

Kitani A, Chua K, Nakamura K, Strober W. Activated self-MHC-reactive T cells have the cytokine phenotype of Th3/T regulatory cell 1 T cells. *J Immunol*. 2000; 165(2):691-702. (Clone-specific: Flow cytometry)

Moran TP, Collier M, McKinnon KP, Davis NL, Johnston RE, Serody JS. A novel viral system for generating antigen-specific T cells. *J Immunol*. 2008; 175(5):3431-3438. (Clone-specific: Flow cytometry)

Sorg RV, Kogler G, Wernet P. Identification of cord blood dendritic cells as an immature CD11c- population. *Blood*. 1999; 93(7):2302-2307. (Clone-specific: Flow cytometry)