

Technical Data Sheet


FITC Mouse Anti-Human TNF

Product Information

Material Number:	562082
Alternate Name:	Tumor necrosis factor alpha; TNF- α ; TNF- α ; TNFSF2; Cachectin
Size:	25 μ g
Concentration:	0.5 mg/ml
Clone:	MAb11
Immunogen:	Recombinant Human TNF
Isotype:	Mouse IgG1, κ
Reactivity:	QC Testing: Human Tested in Development: Rhesus, Cynomolgus, Baboon
Storage Buffer:	Aqueous buffered solution containing \leq 0.09% sodium azide.

Description

The MAb11 monoclonal antibody specifically binds to human tumor necrosis factor (TNF, also known as TNF- α) protein. TNF is an efficient juxtacrine, paracrine and endocrine mediator of inflammatory and immune functions. It regulates the growth and differentiation of a variety of cell types. TNF is cytotoxic for transformed cells when in conjunction with IFN- γ . It is secreted by activated monocytes/macrophages and other cells such as B cells, T cells and fibroblasts. The immunogen used to generate the MAb11 hybridoma was recombinant human TNF. The MAb11 antibody has been reported to crossreact with Rhesus Macaque TNF.

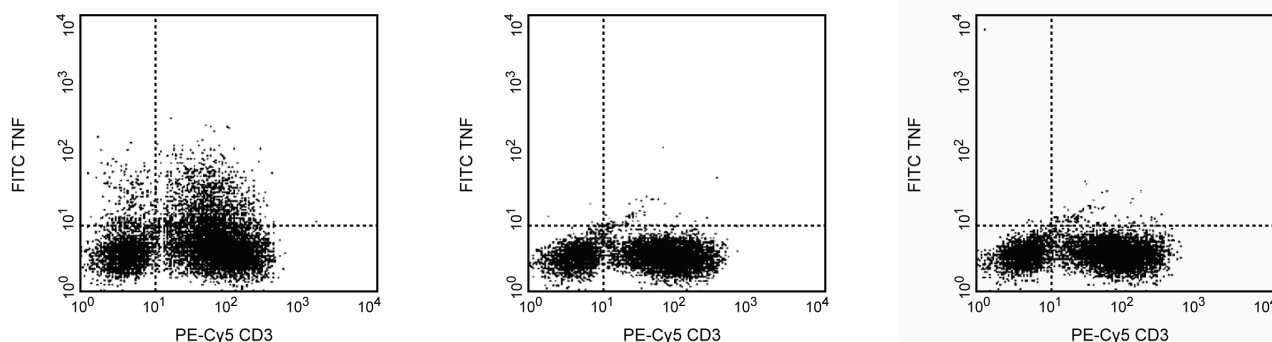

Expression of TNF by stimulated human peripheral blood mononuclear cells (PBMC). Human PBMC were stimulated for 6 hours with PMA (Sigma, Cat. No. P-8139) and calcium ionophore A23187 (Sigma, Cat. No. C-9275) in the presence of GolgiStop™ (aka 2 μ M monensin, Cat. No. 554724). The PBMC were stained with PE-Cy5™5 Mouse Anti-Human CD3 (Cat. No. 555334), fixed, permeabilized, and subsequently stained with 0.25 μ g of FITC Mouse Anti-Human TNF (Cat. No. 562082/554512; left panel). To demonstrate specificity of staining, the binding of FITC-MAb11 was blocked by the preincubation of the conjugated antibody with molar excess of recombinant human TNF (0.5 μ g; Cat. No. 554618; middle panel), and by preincubation of the fixed/permeabilized cells with an excess of Purified Mouse Anti-Human TNF antibody (10 μ g, Cat. No. 554510; right panel) prior to staining with the FITC-MAb11 antibody. The quadrant markers for the bivariate dot plots were set based on the autofluorescence control, and verified with the recombinant cytokine blocking and unlabeled antibody blocking specificity controls.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with FITC under optimum conditions, and unreacted FITC was removed.

Application Notes

Application

Intracellular staining (flow cytometry)	Routinely Tested
---	------------------

Recommended Assay Procedure:

The FITC-conjugated MAb11 antibody can be used for multicolor immunofluorescent staining and flow cytometric analyses to identify and enumerate TNF-producing cells within mixed cell populations (see figure). For optimal immunofluorescent staining with flow cytometric analysis, this anti-cytokine antibody should be titrated (\leq 0.5 μ g mAb/million cells). For specific methodology, please visit the protocols section under "Cytokines (Intracellular Staining)" or "Intracellular Flow" at our website, <http://www.bdbiosciences.com/us/s/resources>.

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

© 2016 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.

562082 Rev. 3



A useful control for demonstrating specificity of staining is either of the following: 1) pre-block the conjugated MAb11 antibody with a molar excess of ligand (e.g., recombinant human TNF; Cat No. 554618) prior to staining, or 2) pre-block the fixed/ permeabilized cells with unlabeled MAb11 antibody (Cat. No. 554510) prior to staining. The staining technique and blocking controls are described in detail by C. Prussin and D. Metcalfe. A suitable mouse IgG1 isotype control for assessing the level of background staining on paraformaldehyde-fixed/saponin-permeabilized human cells is FITC-MOPC-21 (Cat. No. 554679); use at comparable concentrations to antibody of interest.

Suggested Companion Products

<u>Catalog Number</u>	<u>Name</u>	<u>Size</u>	<u>Clone</u>
554618	Recombinant Human TNF	10 µg	(none)
554679	FITC Mouse IgG1, κ Isotype Control	0.1 mg	MOPC-21
554715	BD Cytotfix/Cytoperm Plus Kit (with BD GolgiStop)	250 Tests	(none)
555061	HiCK-1 Human Cytokine Positive Control Cells	1 mL	(none)
554656	Stain Buffer (FBS)	500 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)
554512	FITC Mouse Anti-Human TNF	0.1 mg	MAb11
555334	PE-Cy TM 5 Mouse Anti-Human CD3	100 Tests	UCHT1
554510	Purified Mouse Anti-Human TNF	0.1 mg	MAb11

Product Notices

1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
5. Species testing during development may have been performed with a different format of the same clone. Selected applications have been tested for cross-reactivity.
6. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

- Andersson J, Abrams J, Bjork L, et al. Concomitant in vivo production of 19 different cytokines in human tonsils. *Immunology*. 1994; 83(1):16-24. (Biology)
- Andersson U, Andersson J. Immunolabeling of cytokine-producing cells in tissues and in suspension. In: Fradelizie D, Emelie D, ed. *Cytokine Producing Cells*. Paris: Inserm; 1994:32-49. (Biology)
- Danis VA, Franic GM, Rathjen DA, Brooks PM. Effects of granulocyte-macrophage colony-stimulating factor (GM-CSF), IL-2, interferon-gamma (IFN-gamma), tumour necrosis factor-alpha (TNF-alpha) and IL-6 on the production of immunoreactive IL-1 and TNF-alpha by human monocytes. *Clin Exp Immunol*. 1991; 85(1):143-150. (Clone-specific: ELISA)
- Jason J, Lamed J. Single-cell cytokine profiles in normal humans: comparison of flow cytometric reagents and stimulation protocols. *J Immunol Methods*. 1997; 207(1):13-22. (Biology)
- Litton M, Andersson J, Bjork L, Fehniger T, Ulfgren AK, Andersson U. Cytoplasmic cytokine staining in individual cells. In: Debets and Savelkoul, ed. *Human Cytokine Protocols*. Humana Press; 1996(Biology)
- Prussin C, Metcalfe DD. Detection of intracytoplasmic cytokine using flow cytometry and directly conjugated anti-cytokine antibodies. *J Immunol Methods*. 1995; 188(1):117-128. (Methodology)
- Raqib R, Lindberg AA, Wretling B, Bardhan PK, Andersson U, Andersson J. Persistence of local cytokine production in shigellosis in acute and convalescent stages. *Infect Immun*. 1995; 63(1):289-296. (Biology)
- Rathjen DA, Cowan K, Furphy LJ, Aston R. Antigenic structure of human tumour necrosis factor: recognition of distinct regions of TNF alpha by different tumour cell receptors. *Mol Immunol*. 1991; 28(1-2):79-86. (Clone-specific: ELISA)
- Sander B, Andersson J, Andersson U. Assessment of cytokines by immunofluorescence and the paraformaldehyde-saponin procedure. *Immunol Rev*. 1991; 119:65-93. (Biology)
- Sopper S, Stahl-Hennig C, Demuth M, Johnston IC, Dorries R, ter Meulen V. Lymphocyte subsets and expression of differentiation markers in blood and lymphoid organs of rhesus monkeys. *Cytometry*. 1997; 29(4):351-362. (Biology)
- Verdier F, Aujoulat M, Condevaux F, Descotes J. Determination of lymphocyte subsets and cytokine levels in cynomolgus monkeys. *Toxicology*. 1995; 105(1):81-90. (Biology)