

Technical Data Sheet

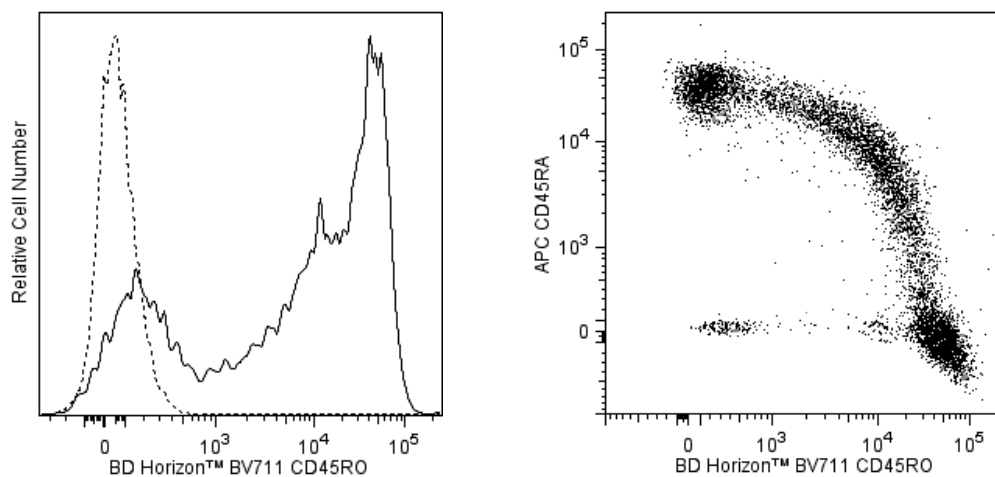
BV711 Mouse Anti-Human CD45RO**Product Information**

Material Number:	563723
Alternate Name:	CD45R; PTPRC; LCA; Leukocyte common antigen; GP180; LY5; T200
Size:	25 Tests
Vol. per Test:	5 µl
Clone:	UCLH1
Immunogen:	Human IL-2-dependent T-cell line
Isotype:	Mouse (BALB/c) IgG2a, κ
Reactivity:	QC Testing: Human
Workshop:	IV N31
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The UCLH1 monoclonal antibody specifically binds to the 180 kDa isoform of CD45 (aka, the Leukocyte Common Antigen). CD45RO is a type I transmembrane glycoprotein that has cytoplasmic protein tyrosine phosphatase activity and functions in signal transduction pathways. This CD45 isoform does not include amino acid sequences encoded by the variable CD45 exons A, B, or C. CD45RO is expressed on most thymocytes, activated T cells, memory T cells, granulocytes and monocytes, but only on a proportion of resting T cells. CD45RO and CD45RA antibodies seem to define complementary, predominantly non-overlapping, populations in resting peripheral T cells, demonstrating heterogeneity within the CD8 and CD4 subpopulations. CD45RO binds to CD22.

The antibody was conjugated to BD Horizon BV711 which is part of the BD Horizon Brilliant™ Violet family of dyes. This dye is a tandem fluorochrome of BD Horizon BV421 with an Ex Max of 405-nm and an acceptor dye with an Em Max at 711-nm. BD Horizon BV711 can be excited by the violet laser and detected in a filter used to detect Cy™5.5 / Alexa Fluor® 700-like dyes (eg, 712/20-nm filter). Due to the excitation and emission characteristics of the acceptor dye, there may be moderate spillover into the Alexa Fluor® 700 and PerCP-Cy5.5 detectors. However, the spillover can be corrected through compensation as with any other dye combination.

**Multicolor flow cytometric analysis of CD45RO expression on human peripheral blood lymphocytes.**

Left Panel: Human whole blood was stained with either BD Horizon™ BV711 Mouse Anti-Human CD45RO antibody (Cat. No. 563722/563723; solid line histogram), or BD Horizon™ BV711 Mouse IgG2a, κ Isotype Control (Cat. No. 563345; dashed line histogram).

Right Panel: Whole blood was stained with BD Horizon™ BV711 Mouse Anti-Human CD45RO and APC Mouse Anti-Human CD45RA (Cat. No. 550855/561884) antibodies.

Erythrocytes were lysed with BD FACS™ Lysing Solution (Cat. No. 349202). The fluorescence histograms and dot plot (CD45RO versus CD45RA expression) were derived from gated events with the forward and side light-scatter characteristics of intact lymphocytes. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.
© 2017 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.



Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ BV711 under optimum conditions, and unconjugated antibody and free BD Horizon™ BV711 were removed.

Application Notes

Application

Flow cytometry

Routinely Tested

Recommended Assay Procedure:

For optimal and reproducible results, BD Horizon Brilliant Stain Buffer should be used anytime two or more BD Horizon Brilliant dyes are used in the same experiment. Fluorescent dye interactions may cause staining artifacts which may affect data interpretation. The BD Horizon Brilliant Stain Buffer was designed to minimize these interactions. More information can be found in the Technical Data Sheet of the BD Horizon Brilliant Stain Buffer (Cat. No. 563794/566349) or the BD Horizon Brilliant Stain Buffer Plus (Cat. No. 566385).

Suggested Companion Products

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 mL	(none)
563345	BV711 Mouse IgG2a, κ Isotype Control	50 µg	G155-178
550855	APC Mouse Anti-Human CD45RA	100 Tests	HI100
561884	APC Mouse Anti-Human CD45RA	25 Tests	HI100
349202	BD FACS™ Lysing Solution	100 mL	(none)
555899	Lysing Buffer	100 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)
563794	Brilliant Stain Buffer	100 Tests	(none)
566349	Brilliant Stain Buffer	1000 Tests	(none)
566385	Brilliant Stain Buffer Plus	1000 Tests	(none)
563722	BV711 Mouse Anti-Human CD45RO	100 Tests	UCHL1

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100-µl experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
5. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
7. BD Horizon Brilliant Violet 711 is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,227,187; 8,455,613; 8,575,303; 8,354,239.
8. BD Horizon Brilliant Stain Buffer is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,575,303; 8,354,239.
9. Cy is a trademark of GE Healthcare.
10. Please refer to www.bdbiosciences.com/pharming/protocols for technical protocols.

References

- Akbar AN, Terry L, Timms A, Beverley PC, Janosy G. Loss of CD45R and gain of UCHL1 reactivity is a feature of primed T cells. *J Immunol.* 1988; 140(7):2171-2178. (Clone-specific: Flow cytometry)
- Knapp W. W. Knapp .. et al., ed. *Leucocyte typing IV : white cell differentiation antigens.* Oxford New York: Oxford University Press; 1989:1-1182. (Clone-specific)
- Norton AJ, Ramsay AD, Smith SH, Beverley PC, Isaacson PG. Monoclonal antibody (UCHL1) that recognises normal and neoplastic T cells in routinely fixed tissues. *J Clin Pathol.* 1986; 39(4):399-405. (Immunogen: Immunohistochemistry)
- Smith SH, Brown MH, Rowe D, Callard RE, Beverley PC. Functional subsets of human helper-inducer cells defined by a new monoclonal antibody, UCHL1. *Immunology.* 1986; 58(1):63-70. (Clone-specific: Flow cytometry, Immunohistochemistry, Immunoprecipitation)
- Streuli M, Morimoto C, Schrieber M, Schlossman SF, Saito H. Characterization of CD45 and CD45R monoclonal antibodies using transfected mouse cell lines that express individual human leukocyte common antigens. *J Immunol.* 1988; 141(11):3910-3914. (Clone-specific: Flow cytometry)
- Zola H. *Leukocyte and stromal cell molecules : the CD markers.* Hoboken, N.J.: Wiley-Liss; 2007(Biology)