

Technical Data Sheet


BV650 Mouse Anti-Human CD62L**Product Information**

Material Number:	563808
Alternate Name:	SELL; L-selectin; LSEL; LAM-1; LECAM-1; LEU8; LNHR; MEL-14; PLNHR; TQ-1
Size:	100 tests
Vol. per Test:	5 µl
Clone:	DREG-56
Immunogen:	Supernatant from PMA-activated Human Peripheral Blood Leukocytes
Isotype:	Mouse IgG1, κ
Reactivity:	QC Testing: Human
Workshop:	V S056
Storage Buffer:	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The DREG-56 monoclonal antibody specifically binds to CD62L. CD62L is a 76-95 kDa glycoprotein that is also referred to as L-selectin or LECAM-1. CD62L is expressed on neutrophils, monocytes, T- and B-lymphocyte subsets and NK cells. The DREG-56 antibody recognizes the same antigen as LAM-1, and specifically inhibits >90% of binding of human lymphocytes to high endothelial venules (HEV) in frozen sections of peripheral, but not mucosal lymphoid tissue. It thus defines L-selectin as a human lymphocyte homing receptor for peripheral lymph node HEV.

The antibody was conjugated to BD Horizon™ BV650 which is part of the BD Horizon™ Brilliant Violet™ family of dyes. This dye is a tandem fluorochrome of BD Horizon™ BV421 with an Ex Max of 405-nm and an acceptor dye with an Em Max at 650-nm. BD Horizon™ BV650 can be excited by the violet laser and detected in a filter used to detect APC-like dyes (eg, 660/20-nm filter). Due to the excitation and emission characteristics of the acceptor dye, there will be spillover into the APC and Alexa Fluor® 700 detectors. However, the spillover can be corrected through compensation as with any other dye combination.

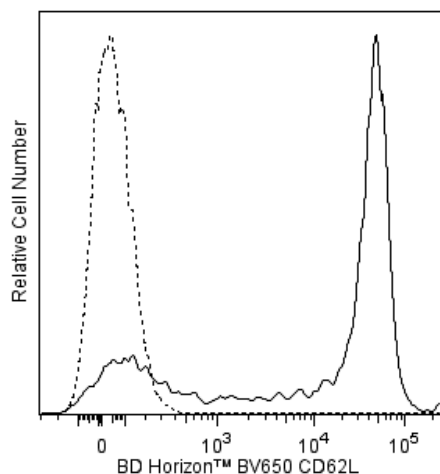

Flow cytometric analysis of CD62L expression on human peripheral blood lymphocytes. Human whole blood was stained with either BD Horizon™ BV650 Mouse Anti-Human CD62L antibody (Cat. No. 563808; solid line histogram) or BD Horizon™ BV650 Mouse IgG1, κ Isotype Control (Cat. No. 563231; dashed line histogram). The erythrocytes were lysed with BD FACS™ Lysing Solution (Cat. No. 349202). The fluorescence histograms were derived from events with the forward and side light-scatter characteristics of intact lymphocytes. Flow cytometric analysis was performed using a BD™ LSR II Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ BV650 under optimum conditions, and unconjugated antibody and free BD Horizon™ BV650 were removed.

Application Notes**Application**

Flow cytometry	Routinely Tested
----------------	------------------

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

Suggested Companion Products

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 ml	(none)
563231	BV650 Mouse IgG1, k Isotype Control	50 µg	X40
349202	BD FACSTM Lysing Solution	100 ml	(none)
555899	Lysing Buffer	100 ml	(none)

Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1×10^6 cells in a 100-µl experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
5. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
6. Brilliant Violet™ 650 is a trademark of Sirigen.
7. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
8. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

Kishimoto TK, Jutila MA, Butcher EC. Identification of a human peripheral lymph node homing receptor: a rapidly down-regulated adhesion molecule. *Proc Natl Acad Sci U S A*. 1990; 87(6):2244-2248. (Immunogen: Flow cytometry, Inhibition)

Kishimoto TK, Warnock RA, Jutila MA, et al. Antibodies against human neutrophil LECAM-1 (LAM-1/Leu-8/DREG-56 antigen) and endothelial cell ELAM-1 inhibit a common CD18-independent adhesion pathway in vitro. *Blood*. 1991; 78(3):805-811. (Clone-specific: Inhibition)

Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. New York: Oxford University Press; 1995. (Clone-specific: Flow cytometry, Immunocytochemistry (cytospins), Immunohistochemistry)

Zola H, Swart B, Nicholson I, Voss E. *Leukocyte and Stromal Cell Molecules. The CD Markers*. Hoboken, New Jersey: John Wiley & Sons, Inc.; 2007:1-581. (Biology)

BD Biosciences

bdbiosciences.com

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	800.979.9408	32.53.720.550	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

