

## Technical Data Sheet

## BUV737 Rat Anti-Mouse CD8a


## Product Information

<b>Material Number:</b>	564297
<b>Alternate Name:</b>	CD8a; CD8 alpha chain; Ly-2; Lyt2; Lyt-2; Ly-35; Ly-B
<b>Size:</b>	50 µg
<b>Concentration:</b>	0.2 mg/ml
<b>Clone:</b>	53-6.7
<b>Immunogen:</b>	Mouse Spleen Cells or Thymocyte Membranes
<b>Isotype:</b>	Rat (LOU) IgG2a, κ
<b>Reactivity:</b>	QC Testing: Mouse
<b>Storage Buffer:</b>	Aqueous buffered solution containing ≤0.09% sodium azide.

## Description

The 53-6.7 monoclonal antibody specifically binds to the 38 kDa  $\alpha$  and 34 kDa  $\alpha'$  chains of the CD8 differentiation antigen (Ly-2 or Lyt-2) of all mouse strains tested. The CD8  $\alpha$  and  $\alpha'$  chains (CD8a) form heterodimers with the CD8  $\beta$  chain (CD8b, Ly-3, or Lyt-3) on the surface of most thymocytes. A subpopulation of mature T lymphocytes (i.e., MHC class I-restricted T cells, including most T suppressor/cytotoxic cells) expresses almost exclusively the CD8  $\alpha\beta$  heterodimer. Subsets of  $\gamma\delta$  TCR-bearing T cells, intestinal intraepithelial lymphocytes, and dendritic cells express CD8a without CD8b. It has been suggested that the expression of the CD8a/CD8b heterodimer is restricted to T lymphocytes which matured in the thymus or in an extrathymic environment that had been influenced by thymus-initiated neuroendocrine signals. CD8 is an antigen coreceptor on the T-cell surface which interacts with MHC class I molecules on antigen-presenting cells or epithelial cells. It participates in T-cell activation through its association with the T-cell receptor complex and protein tyrosine kinase lck (p56 [lck]). The CD8  $\alpha$  and  $\alpha'$  chains arise from alternatively spliced messengers of a single *CD8a* gene. The longer  $\alpha$  form associates with p56 [lck] via a CXCP motif in its cytoplasmic domain, which it shares with CD4, but not with CD8b. The truncated  $\alpha'$  chain is unable to associate with p56 [lck], and it may function to attenuate the CD8-mediated costimulatory signal during intrathymic T-cell maturation. In vivo and in vitro treatment with 53-6.7 mAb has reportedly been effective at depleting CD8<sup>+</sup> peripheral T lymphocytes. The 53-6.7 antibody has also been reported to cross-react with CD8  $\alpha$ - and  $\alpha'$ -like polypeptides on subsets of thymic and peripheral lymphocytes in the Egyptian toad, *Bufo regularis*.

The antibody was conjugated to BD Horizon BUV737 which is part of the BD Horizon Brilliant™ Ultraviolet family of dyes. This dye is a tandem fluorochrome with an Ex Max near 350 nm and an Em Max near 737 nm. BD Horizon Brilliant BUV737 can be excited by the ultraviolet laser (355 nm) and detected with a 740/35 nm filter. Due to the excitation of the acceptor dye by the red laser line, there may be significant spillover into red laser detectors with filters in the 700-720 nm range.

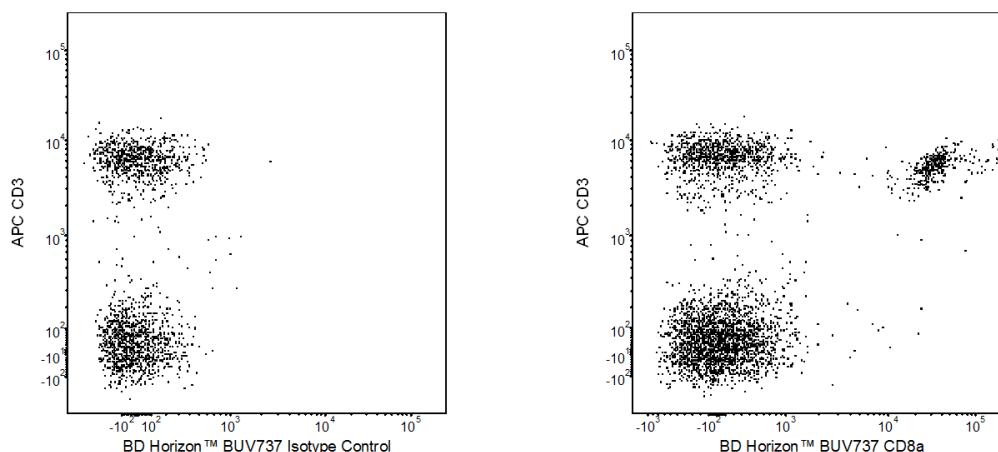

**Two-color flow cytometric analysis of CD8a expression on mouse splenocytes.** Mouse splenic leucocytes were preincubated with Purified Rat Anti-Mouse CD16/CD32 antibody (Mouse BD Fc Block™) (Cat. No. 553141/553142). The cells were then stained with APC Hamster Anti-Mouse CD3e (Cat. No. 561826/553066) and either BD Horizon™ BUV737 Rat IgG2a, κ Isotype Control (Cat. No. 564294; Left Panel) or BD Horizon BUV737 Rat Anti-Mouse CD8a antibody (Cat. No. 564297; Right Panel) antibodies. The two-color fluorescence dot plot shows the correlated expression patterns of CD8a versus CD3e for gated events with the forward and side light-scatter characteristic of viable splenic leucocytes. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

## BD Biosciences

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.  
© 2017 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.



## Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon BUV737 under optimum conditions, and unconjugated antibody and free BD Horizon BUV737 were removed.

## Application Notes

### Application

Flow cytometry

Routinely Tested

### Recommended Assay Procedure:

BD™ CompBeads can be used as surrogates to assess fluorescence spillover (Compensation). When fluorochrome conjugated antibodies are bound to BD CompBeads, they have spectral properties very similar to cells. However, for some fluorochromes there can be small differences in spectral emissions compared to cells, resulting in spillover values that differ when compared to biological controls. It is strongly recommended that when using a reagent for the first time, users compare the spillover on cells and BD CompBead to ensure that BD CompBeads are appropriate for your specific cellular application.

For optimal and reproducible results, BD Horizon Brilliant Stain Buffer should be used anytime two or more BD Horizon Brilliant dyes are used in the same experiment. Fluorescent dye interactions may cause staining artifacts which may affect data interpretation. The BD Horizon Brilliant Stain Buffer was designed to minimize these interactions. More information can be found in the Technical Data Sheet of the BD Horizon Brilliant Stain Buffer (Cat. No. 563794/566349) or the BD Horizon Brilliant Stain Buffer Plus (Cat. No. 566385).

Note: When using high concentrations of antibody, background binding of this dye to erythroid cell subsets (mature erythrocytes and precursors) has been observed. For researchers studying these cell populations, or in cases where light scatter gating does not adequately exclude these cells from the analysis, this background may be an important factor to consider when selecting reagents for panel(s).

### Suggested Companion Products

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)
564294	BUV737 Rat IgG2a, κ Isotype Control	50 µg	R35-95
555899	Lysing Buffer	100 mL	(none)
561826	APC Hamster Anti-Mouse CD3e	25 µg	145-2C11
553066	APC Hamster Anti-Mouse CD3e	0.1 mg	145-2C11
553141	Purified Rat Anti-Mouse CD16/CD32 (Mouse BD Fc Block™)	0.1 mg	2.4G2
553142	Purified Rat Anti-Mouse CD16/CD32 (Mouse BD Fc Block™)	0.5 mg	2.4G2
349202	BD FACSTM Lysing Solution	100 mL	(none)
563794	Brilliant Stain Buffer	100 Tests	(none)
566349	Brilliant Stain Buffer	1000 Tests	(none)
566385	Brilliant Stain Buffer Plus	1000 Tests	(none)

### Product Notices

1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
5. BD Horizon Brilliant Ultraviolet 737 is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,227,187; 8,575,303; 8,354,239.
6. BD Horizon Brilliant Stain Buffer is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,575,303; 8,354,239.
7. Please refer to <http://regdocs.bd.com> to access safety data sheets (SDS).
8. Please refer to [www.bdbiosciences.com/pharming/en/protocols](http://www.bdbiosciences.com/pharming/en/protocols) for technical protocols.

### References

Bierer BE, Sleckman BP, Ratnofsky SE, Burakoff SJ. The biologic roles of CD2, CD4, and CD8 in T-cell activation. *Annu Rev Immunol.* 1989; 7:579-599. (Biology)

Fujiura Y, Kawaguchi M, Kondo Y, et al. Development of CD8 alpha alpha+ intestinal intraepithelial T cells in beta 2-microglobulin- and/or TAP1-deficient mice. *J Immunol.* 1996; 156(8):2710-2715. (Clone-specific; Flow cytometry)

Hathcock KS. T cell depletion by cytotoxic elimination. *Curr Protoc Immunol.* 1991; 1:3.4.1-3.4.3. (Clone-specific; Cell separation, Depletion, Flow cytometry)

Kruisbeek AM, Shevach EM. Proliferative assays for T cell function. *Curr Protoc Immunol.* 2004; 3:3.12.1-3.12.14. (Clone-specific; Depletion, In vivo exacerbation)

Ledbetter JA, Herzenberg LA. Xenogeneic monoclonal antibodies to mouse lymphoid differentiation antigens. *Immunol Rev.* 1979; 47:63-90. (Immunogen; Flow cytometry, Immunofluorescence, Immunoprecipitation)

Ledbetter JA, Rouse RV, Micklem HS, Herzenberg LA. T cell subsets defined by expression of Lyt-1,2,3 and Thy-1 antigens. Two-parameter immunofluorescence and cytotoxicity analysis with monoclonal antibodies modifies current views. *J Exp Med.* 1980; 152(2):280-295. (Immunogen; Flow cytometry)

Ledbetter JA, Seaman WE, Tsu TT, Herzenberg LA. Lyt-2 and lyt-3 antigens are on two different polypeptide subunits linked by disulfide bonds. Relationship of subunits to T cell cytolytic activity. *J Exp Med.* 1981; 153(6):1503-1516. (Clone-specific: Blocking, Flow cytometry, Immunoprecipitation, Inhibition)

LeFrancois L. Extrathymic differentiation of intraepithelial lymphocytes: generation of a separate and unequal T-cell repertoire. *Immunol Today.* 1991; 12(12):436-438. (Biology)

Leishman AJ, Naidenko OV, Attinger A, et al. T cell responses modulated through interaction between CD8alphaalpha and the nonclassical MHC class I molecule, TL. *Science.* 2001; 294(5548):1848-1849. (Clone-specific: Blocking, Flow cytometry)

MacDonald HR, Schreyer M, Howe RC, Bron C. Selective expression of CD8 alpha (Lyt-2) subunit on activated thymic gamma/delta cells. *Eur J Immunol.* 1990; 20(4):927-930. (Biology)

Nakayama K, Nakayama K, Negishi I, et al. Requirement for CD8 beta chain in positive selection of CD8-lineage T cells. *Science.* 1994; 263(5150):1131-1133. (Biology)

Sydora BC, Brossay L, Hagenbaugh A, Kronenberg M, Cheroutre H. TAP-independent selection of CD8+ intestinal intraepithelial lymphocytes. *J Immunol.* 1996; 156(11):4209-4216. (Clone-specific: Flow cytometry)

Takahashi K, Nakata M, Tanaka T, et al. CD4 and CD8 regulate interleukin 2 responses of T cells. *Proc Natl Acad Sci U S A.* 1992; 89(12):5557-5561. (Clone-specific: Immunoprecipitation, Inhibition)

van Ewijk W, van Soest PL, van den Engh GJ. Fluorescence analysis and anatomic distribution of mouse T lymphocyte subsets defined by monoclonal antibodies to the antigens Thy-1, Lyt-1, Lyt-2, and T-200. *J Immunol.* 1981; 127(6):2594-2604. (Clone-specific: Flow cytometry, Immunofluorescence, Immunohistochemistry)

Vremec D, Zorbas M, Scollay R, et al. The surface phenotype of dendritic cells purified from mouse thymus and spleen: investigation of the CD8 expression by a subpopulation of dendritic cells. *J Exp Med.* 1992; 176(1):47-58. (Clone-specific: Flow cytometry)

Wang J, Klein JR. Thymus-neuroendocrine interactions in extrathymic T cell development. *Science.* 1994; 265(5180):1860-1862. (Biology)

Wu L, Vremec D, Ardavin C, et al. Mouse thymus dendritic cells: kinetics of development and changes in surface markers during maturation. *Eur J Immunol.* 1995; 25(2):418-425. (Biology)

Zamoyska R, Derham P, Gorman SD, et al. Inability of CD8 alpha' polypeptides to associate with p56lck correlates with impaired function in vitro and lack of expression in vivo. *Nature.* 1989; 342(6247):278-281. (Biology)