

## Technical Data Sheet


**BUV737 Mouse Anti-Human CD20****Product Information**

<b>Material Number:</b>	<b>564431</b>
<b>Alternate Name:</b>	MS4A1; B1; Bp35; LEU-16; S7
<b>Size:</b>	25 Tests
<b>Vol. per Test:</b>	5 µl
<b>Clone:</b>	2H7
<b>Immunogen:</b>	Human 6.16c1.3 B cell line
<b>Isotype:</b>	Mouse (C57BL/6) IgG2b, κ
<b>Reactivity:</b>	QC Testing: Human Tested in Development: Rhesus, Cynomolgus, Baboon
<b>Workshop:</b>	II B22; III B739, NL382; IV B201
<b>Storage Buffer:</b>	Aqueous buffered solution containing ≤0.09% sodium azide.

**Description**

The 2H7 monoclonal antibody specifically binds to CD20, encoded by the *MS4A1* (Membrane-spanning 4-domains, subfamily A, member 1) gene. CD20 is a 33-37 kDa, unglycosylated four-transmembrane phosphoprotein. CD20 is expressed on pre-B-cells, resting and activated B cells, and follicular dendritic cells, but not plasma cells. Low level CD20 expression is observed on a small subset of normal circulating T lymphocytes. The CD20 molecule is involved in the regulation of B-cell activation.

The antibody was conjugated to BD Horizon BUV737 which is part of the BD Horizon Brilliant™ Ultraviolet family of dyes. This dye is a tandem fluorochrome with an Ex Max near 350 nm and an Em Max near 737 nm. BD Horizon Brilliant BUV737 can be excited by the ultraviolet laser (355 nm) and detected with a 740/35 nm filter. Due to the excitation of the acceptor dye by the red laser line, there may be significant spillover into red laser detectors with filters in the 700-720 nm range.

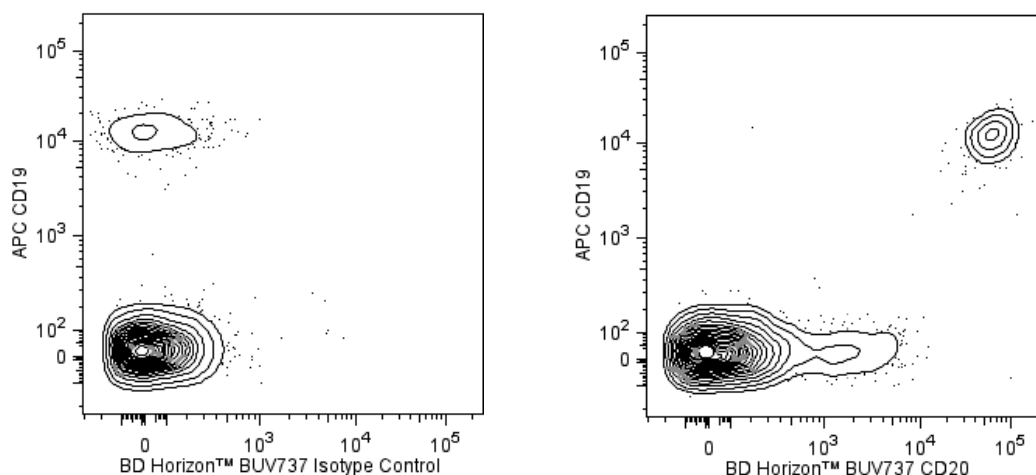

**Two-color flow cytometric analysis of CD20 expression on human peripheral blood lymphocytes.** Whole blood was stained with APC Mouse Anti-Human CD19 antibody (Cat. No. 555415/561742) and either BD Horizon™ BUV737 Mouse IgG2b, κ Isotype Control (Cat. No. 564429; Left Panel) or BD Horizon BUV737 Mouse Anti-Human CD20 antibody (Cat. No. 564431/564432; Right Panel). Erythrocytes were lysed with BD FACS Lysing Solution (Cat. No. 349202). Two color flow cytometric contour plots showing the correlated expression of CD20 (or Ig Isotype Control staining) versus CD19 were derived from gated events with the forward light-scatter characteristics of intact lymphocytes. Flow cytometry was performed using a BD™ LSR II Flow Cytometer System.

**Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon BUV737 under optimum conditions, and unconjugated antibody and free BD Horizon BUV737 were removed.

**BD Biosciences**

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

*Conditions:* The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.  
© 2017 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.



## Application Notes

### Application

Flow cytometry

Routinely Tested

### Recommended Assay Procedure:

BD™ CompBeads can be used as surrogates to assess fluorescence spillover (Compensation). When fluorochrome conjugated antibodies are bound to BD CompBeads, they have spectral properties very similar to cells. However, for some fluorochromes there can be small differences in spectral emissions compared to cells, resulting in spillover values that differ when compared to biological controls. It is strongly recommended that when using a reagent for the first time, users compare the spillover on cells and BD CompBead to ensure that BD CompBeads are appropriate for your specific cellular application.

For optimal and reproducible results, BD Horizon Brilliant Stain Buffer should be used anytime two or more BD Horizon Brilliant dyes are used in the same experiment. Fluorescent dye interactions may cause staining artifacts which may affect data interpretation. The BD Horizon Brilliant Stain Buffer was designed to minimize these interactions. More information can be found in the Technical Data Sheet of the BD Horizon Brilliant Stain Buffer (Cat. No. 563794/566349) or the BD Horizon Brilliant Stain Buffer Plus (Cat. No. 566385).

Note: When using high concentrations of antibody, background binding of this dye to erythroid cell subsets (mature erythrocytes and precursors) has been observed. For researchers studying these cell populations, or in cases where light scatter gating does not adequately exclude these cells from the analysis, this background may be an important factor to consider when selecting reagents for panel(s).

### Suggested Companion Products

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)
564429	BUV737 Mouse IgG2b, κ Isotype Control	50 µg	27-35
349202	BD FACSTM Lysing Solution	100 mL	(none)
555899	Lysing Buffer	100 mL	(none)
564432	BUV737 Mouse Anti-Human CD20	100 Tests	2H7
555415	APC Mouse Anti-Human CD19	100 Tests	HIB19
561742	APC Mouse Anti-Human CD19	25 Tests	HIB19
563794	Brilliant Stain Buffer	100 Tests	(none)
566349	Brilliant Stain Buffer	1000 Tests	(none)
566385	Brilliant Stain Buffer Plus	1000 Tests	(none)

### Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100-µl experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. Species testing during development may have been performed with a different format of the same clone. Selected applications have been tested for cross-reactivity.
5. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
6. Please refer to <http://regdocs.bd.com> to access safety data sheets (SDS).
7. BD Horizon Brilliant Stain Buffer is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,575,303; 8,354,239.
8. BD Horizon Brilliant Ultraviolet 737 is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,227,187; 8,575,303; 8,354,239.
9. Please refer to [www.bdbiosciences.com/pharmingen/protocols](http://www.bdbiosciences.com/pharmingen/protocols) for technical protocols.

### References

Clark EA, Yokochi T. Human B cell and B cell blast-associated surface molecules defined with monoclonal antibodies. In: Bernard A, Boumsell L, Dausset J, Milstein C, Schlossman SF, ed. *Leukocyte Typing*. Berlin: Springer-Verlag; 1984:339-346. (Clone-specific: Blocking, Flow cytometry, Immunoprecipitation)

Hultin LE, Hausner MA, Hultin PM, Giorgi JV. CD20 (pan-B cell) antigen is expressed at a low level on a subpopulation of human T lymphocytes. *Cytometry*. 1993; 14(2):193-204. (Biology)

Knapp W. W. Knapp .. et al., ed. *Leukocyte typing IV : white cell differentiation antigens*. Oxford New York: Oxford University Press; 1989:1-1182. (Biology)

Ledbetter JA, Clark EA. Surface phenotype and function of tonsillar germinal center and mantle zone B cell subsets. *Hum Immunol*. 1986; 15:30-43. (Immunogen: Blocking, Flow cytometry)

Loken MR, Shah VO, Dattilio KL, Civin CI. Flow cytometric analysis of human bone marrow. II. Normal B lymphocyte development. *Blood*. 1987; 70(5):1316-1324. (Biology)

Reinherz EL, Ellis L, Reinherz .. et al., ed. *Leukocyte typing II*. New York: Springer-Verlag; 1986:1-560. (Clone-specific: Flow cytometry, Immunoprecipitation)

Schlossman SF, Stuart F, Schlossman .. et al., ed. *Leukocyte typing V : white cell differentiation antigens : proceedings of the fifth international workshop and conference held in Boston, USA, 3-7 November, 1993*. Oxford: Oxford University Press; 1995(Biology)

Tedder TF, Penta A. Structure of the CD20 antigen and gene of human and mouse B cells: use of transfected cell lines to examine the Workshop panel of antibodies. In: Knapp W. W. Knapp .. et al., ed. *Leukocyte typing IV : white cell differentiation antigens*. Oxford New York: Oxford University Press; 1989:48-50. (Clone-specific: Blocking, Flow cytometry)