

## Technical Data Sheet

## BUV737 Mouse Anti-Human CD56


## Product Information

<b>Material Number:</b>	564447
<b>Alternate Name:</b>	NCAM1; NCAM-1; NCAM; Leu-19; Neural cell adhesion molecule 1; NKH1; MSK39
<b>Size:</b>	100 Tests
<b>Vol. per Test:</b>	5 µl
<b>Clone:</b>	NCAM16.2 (also known as NCAM 16)
<b>Immunogen:</b>	Immunoaffinity-enriched adult human brain NCAM
<b>Isotype:</b>	Mouse (BALB/c) IgG2b, κ
<b>Reactivity:</b>	QC Testing: Human
<b>Workshop:</b>	V NK60
<b>Storage Buffer:</b>	Aqueous buffered solution containing ≤0.09% sodium azide.

## Description

The NCAM16.2 monoclonal antibody specifically binds to human CD56. It recognizes an extracellular immunoglobulin-like domain common to 120, 140, and 180 kDa forms of CD56, also known as the neural cell adhesion molecule (NCAM), NKH1 or MSK39. The CD56 antigen is expressed on approximately 10% to 25% of peripheral blood lymphocytes. It is present on essentially all resting and activated CD16+ natural killer (NK) lymphocytes and approximately 5% of CD3+ peripheral blood lymphocytes. CD3+ CD56+ T lymphocytes comprise a unique subset of cytotoxic T lymphocytes that mediates non-major histocompatibility complex (MHC)-restricted cytotoxicity. CD56 antigen density on NK lymphocytes increases upon cellular activation. The CD56 antigen is involved in neuronal homotypic cell adhesion and cell differentiation during embryogenesis. CD16+ CD56+ NK cells demonstrate reciprocal transfer of an activation state with dendritic cells.

The antibody was conjugated to BD Horizon BUV737 which is part of the BD Horizon Brilliant™ Ultraviolet family of dyes. This dye is a tandem fluorochrome with an Ex Max near 350 nm and an Em Max near 737 nm. BD Horizon Brilliant BUV737 can be excited by the ultraviolet laser (355 nm) and detected with a 740/35 nm filter. Due to the excitation of the acceptor dye by the red laser line, there may be significant spillover into red laser detectors with filters in the 700-720 nm range.

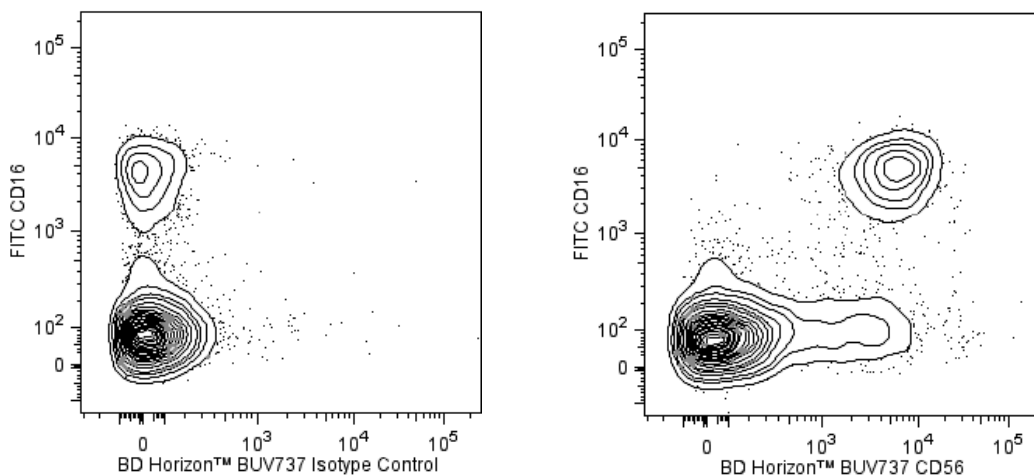

**Two-color flow cytometric analysis of CD56 expression on human peripheral blood lymphocytes.** Human whole blood was stained with FITC Mouse Anti-Human CD16 antibody (Cat. No. 555406/556618/560996) and either BD Horizon™ BUV737 Mouse IgG2b, κ Isotype Control (Cat. No. 564429; Left Panel) or BD Horizon BUV737 Mouse Anti-Human CD56 antibody (Cat. No. 564447/564448; Right Panel). The erythrocytes were lysed with BD FACS™ Lysing Solution (Cat. No. 349202). Two-color flow cytometric dot plots show the correlated expression patterns of CD56 (or Ig Isotype control staining) versus CD16 for gated events with the forward and side light-scatter characteristics of intact lymphocytes. Flow cytometric analysis was performed using a BD LSRFortessa™ Cell Analyzer System.

## BD Biosciences

bdbiosciences.com

United States 877.232.8995   Canada 866.979.9408   Europe 32.2.400.98.95   Japan 0120.8555.90   Asia Pacific 65.6861.0633   Latin America/Caribbean 55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.  
© 2017 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.



## Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon BUV737 under optimum conditions, and unconjugated antibody and free BD Horizon BUV737 were removed.

## Application Notes

### Application

Flow cytometry

Routinely Tested

### Recommended Assay Procedure:

BD™ CompBeads can be used as surrogates to assess fluorescence spillover (Compensation). When fluorochrome conjugated antibodies are bound to BD CompBeads, they have spectral properties very similar to cells. However, for some fluorochromes there can be small differences in spectral emissions compared to cells, resulting in spillover values that differ when compared to biological controls. It is strongly recommended that when using a reagent for the first time, users compare the spillover on cells and BD CompBead to ensure that BD CompBeads are appropriate for your specific cellular application.

For optimal and reproducible results, BD Horizon Brilliant Stain Buffer should be used anytime two or more BD Horizon Brilliant dyes are used in the same experiment. Fluorescent dye interactions may cause staining artifacts which may affect data interpretation. The BD Horizon Brilliant Stain Buffer was designed to minimize these interactions. More information can be found in the Technical Data Sheet of the BD Horizon Brilliant Stain Buffer (Cat. No. 563794/566349) or the BD Horizon Brilliant Stain Buffer Plus (Cat. No. 566385).

Note: When using high concentrations of antibody, background binding of this dye to erythroid cell subsets (mature erythrocytes and precursors) has been observed. For researchers studying these cell populations, or in cases where light scatter gating does not adequately exclude these cells from the analysis, this background may be an important factor to consider when selecting reagents for panel(s).

### Suggested Companion Products

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)
564429	BUV737 Mouse IgG2b, κ Isotype Control	50 µg	27-35
555406	FITC Mouse Anti-Human CD16	100 Tests	3G8
556618	FITC Mouse Anti-Human CD16	50 Tests	3G8
560996	FITC Mouse Anti-Human CD16	25 Tests	3G8
349202	BD FACS™ Lysing Solution	100 mL	(none)
555899	Lysing Buffer	100 mL	(none)
563794	Brilliant Stain Buffer	100 Tests	(none)
564448	BUV737 Mouse Anti-Human CD56	25 Tests	NCAM16.2
566349	Brilliant Stain Buffer	1000 Tests	(none)
566385	Brilliant Stain Buffer Plus	1000 Tests	(none)

### Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100-µl experimental sample (a test).
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
5. Please refer to <http://regdocs.bd.com> to access safety data sheets (SDS).
6. BD Horizon Brilliant Stain Buffer is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,575,303; 8,354,239.
7. BD Horizon Brilliant Ultraviolet 737 is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,227,187; 8,575,303; 8,354,239.
8. Please refer to [www.bdbiosciences.com/pharming/protocols](http://www.bdbiosciences.com/pharming/protocols) for technical protocols.

### References

Bennett IM, Zatzepina O, Zamai L, Azzoni L, Mikheeva T, Perussia B. Definition of a natural killer NKR-P1A+/CD56-/CD16- functionally immature human NK cell subset that differentiates in vitro in the presence of interleukin 12. *J Exp Med*. 1996; 184(5):1845-1856. (Biology)

Campbell JJ, Qin S, Unutmaz D, et al. Unique subpopulations of CD56+ NK and NK-T peripheral blood lymphocytes identified by chemokine receptor expression repertoire. *J Immunol*. 2001; 166(11):6477-6482. (Biology)

Cooper MA, Fehniger TA, Caligiuri MA. The biology of human natural killer-cell subsets. *Trends Immunol*. 2001; 22(11):633-640. (Biology)

Cunningham BA, Hemperly JJ, Murray BA, Prediger EA, Brackenbury R, Edelman GM. Neural cell adhesion molecule: structure, immunoglobulin-like domains, cell surface modulation, and alternative RNA splicing. *Science*. 1987; 236(4803):799-806. (Biology)

Edelman GM. Cell adhesion molecules in the regulation of animal form and tissue pattern. *Annu Rev Cell Biol*. 1986; 2:81-116. (Biology)

Galandrini R, Tassi I, Mattia G, et al. SH2-containing inositol phosphatase (SHIP-1) transiently translocates to raft domains and modulates CD16-mediated cytotoxicity in human NK cells. *Blood*. 2001; 100(13):4581-4589. (Biology)

Gerosa F, Baldani-Guerra B, Nisii C, Marchesini V, Carra G, Trinchieri G. Reciprocal activating interaction between natural killer cells and dendritic cells. *J Exp Med.* 2002; 195(3):327-333. (Biology)

Lanier LL, Chang C, Azuma M, Ruitenberg JJ, Hemperly JJ, Phillips JH. Molecular and functional analysis of human natural killer cell-associated neural cell adhesion molecule (N-CAM/CD56). *J Immunol.* 1991; 146(12):4421-4426. (Immunogen: ELISA, Flow cytometry)

Lanier LL, Le AM, Civin CI, Loken MR, Phillips JH. The relationship of CD16 (Leu-11) and Leu-19 (NKH-1) antigen expression on human peripheral blood NK cells and cytotoxic T lymphocytes. *J Immunol.* 1986; 136(12):4480-4486. (Biology)

Lanier LL, Testi R, Bindl J, Phillips JH. Identity of Leu-19 (CD56) leukocyte differentiation antigen and neural cell adhesion molecule. *J Exp Med.* 1989; 169(6):2233-2238. (Biology)

Nitta T, Yagita H, Sato K, Okumura K. Involvement of CD56 (NKH-1/Leu-19 antigen) as an adhesion molecule in natural killer-target cell interaction. *J Exp Med.* 1989; 170(5):1757-1761. (Biology)

Ritz J, Trinchieri G, Lanier LL. NK-cell Antigens: Section Report. In: Schlossman SF, Stuart F, Schlossman .. et al., ed. *Leucocyte typing V : white cell differentiation antigens : proceedings of the fifth international workshop and conference held in Boston, USA, 3-7 November, 1993.* Oxford: Oxford University Press; 1995:1367-1372. (Clone-specific: Flow cytometry)

Phillips JH, Lanier LL. Dissection of the lymphokine-activated killer phenomenon: relative contribution of peripheral blood natural killer cells and T lymphocytes to cytotoxicity. *J Exp Med.* 1986; 164(3):814-825. (Biology)

Schubert W, Zimmermann K, Cramer M, Starzinski-Powitz A. Lymphocyte antigen Leu-19 as a molecular marker of regeneration in human skeletal muscle. *Proc Natl Acad Sci U S A.* 1989; 86(1):307-311. (Biology)