

## Technical Data Sheet


## PE Rat Anti-Human CX3CR1

## Product Information

<b>Material Number:</b>	<b>565798</b>
<b>Alternate Name:</b>	CCRL1; CMKBRL1; CMKDR1; GPR13; GPRV28; V28; Fractalkine Receptor
<b>Size:</b>	25 Tests
<b>Vol. per Test:</b>	5 µl
<b>Clone:</b>	2A9-1
<b>Immunogen:</b>	Human CX3CR1 Recombinant Protein
<b>Isotype:</b>	Rat IgG2b, κ
<b>Reactivity:</b>	QC Testing: Human
<b>Storage Buffer:</b>	Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

## Description

The 2A9-1 monoclonal antibody specifically binds to human CX3CR1, which is also known as chemokine (C-C) receptor-like 1 (CCRL1), Beta chemokine receptor-like 1 (CMK-BRL-1), G protein-coupled receptor 13 (GPR13), or GPRV28 (V28). CX3CR1 is a seven transmembrane G protein coupled receptor that is expressed by NK cells, T cells, and monocytes. The cellular expression of CX3CR1 is correlated with high levels of intracellular perforin and granzyme B. CX3CR1 serves as a receptor for fractalkine (CX3CL1). Fractalkine is a transmembrane chemokine of the CX3C family that is expressed on activated endothelial cells, neurons, and astrocytes. Interaction of CX3CR1 with fractalkine initiates cellular adhesive and chemotactic responses.

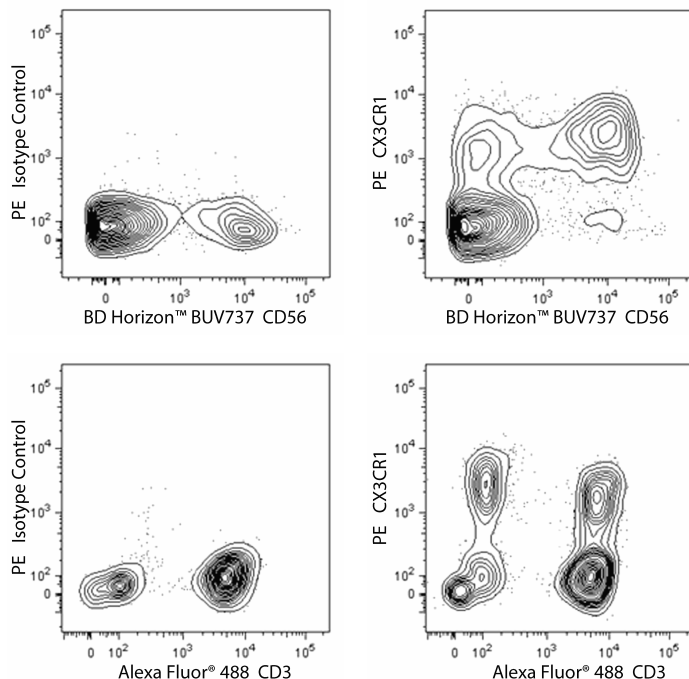

**Multicolor flow cytometric analysis of CX3CR1 expression on human peripheral blood cells.** Whole blood was treated with BD Pharm Lyse™ Lysing Buffer (Cat. No. 555899) to lyse erythrocytes. After washing, the leucocytes were stained with BD Horizon™ BUV737 Mouse Anti-Human CD56 antibody (Cat. No. 564447; Top Plots), Alexa Fluor® 488 Mouse Anti-Human CD3 antibody (Cat. No. 557694; Bottom Plots), and either PE Rat IgG2b, κ Isotype Control (Cat. No. 553989; Left Plots) or PE Rat Anti-Human CX3CR1 antibody (Cat. No. 565796/565798; Right Plots). Two-color flow cytometric contour plots showing the correlated expression of CD56 or CD3 versus CX3CR1 (or Ig Isotype control staining), were derived from gated events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometric analysis was performed using a BD LSRFortessa™ X-20 Flow Cytometer System.

## Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with R-PE under optimum conditions, and unconjugated antibody and free PE were removed.

## Application Notes

## Application

Flow cytometry	Routinely Tested
----------------	------------------

## BD Biosciences

bdbiosciences.com

<b>United States</b>	<b>Canada</b>	<b>Europe</b>	<b>Japan</b>	<b>Asia Pacific</b>	<b>Latin America/Caribbean</b>
877.232.8995	866.979.9408	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2015 BD

565798 Rev. 1



## Suggested Companion Products

Catalog Number	Name	Size	Clone
553989	PE Rat IgG2b, $\kappa$ Isotype Control	0.1 mg	A95-1
557694	Alexa Fluor® 488 Mouse Anti-Human CD3	100 Tests	UCHT1
564447	BUV737 Mouse Anti-Human CD56	100 Tests	NCAM16.2
565796	PE Rat Anti-Human CX3CR1	100 Tests	2A9-1
555899	Lysing Buffer	100 mL	(none)
349202	BD FACS™ Lysing Solution	100 mL	(none)
554656	Stain Buffer (FBS)	500 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)

## Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100- $\mu$ l experimental sample (a test).
2. Please refer to [www.bdbiosciences.com/pharming/en/protocols](http://www.bdbiosciences.com/pharming/en/protocols) for technical protocols.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
5. Alexa Fluor® is a registered trademark of Molecular Probes, Inc., Eugene, OR.
6. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
7. An isotype control should be used at the same concentration as the antibody of interest.

## References

Kobayashi T, Okamoto S, Iwakami Y, et al. Exclusive increase of CX3CR1+CD28-CD4+ T cells in inflammatory bowel disease and their recruitment as intraepithelial lymphocytes. *Inflamm Bowel Dis.* 2007; 13(7):837-46. (Biology)

Kondo Y, Kimura O, Tanaka Y, et al. Differential Expression of CX3CL1 in Hepatitis B Virus-Replicating Hepatoma Cells Can Affect the Migration Activity of CX3CR1+ Immune Cells. *J Virol.* 2015; 89(14):7016-27. (Biology)

Nanki T, Imai T, Nagasaka K, et al. Migration of CX3CR1-positive T cells producing type 1 cytokines and cytotoxic molecules into the synovium of patients with rheumatoid arthritis. *Arthritis Rheum.* 2002; 46(11):2878-83. (Clone-specific: Flow cytometry)

Nishimura M, Umehara H, Nakayama T, et al. Dual functions of fractalkine/CX3C ligand 1 in trafficking of perforin+/granzyme B+ cytotoxic effector lymphocytes that are defined by CX3CR1 expression. *J Immunol.* 2002; 168(12):6173-80. (Immunogen: Flow cytometry)

Siwetz M, Sundl M, Kolb D, et al. Placental fractalkine mediates adhesion of THP-1 monocytes to villous trophoblast. *Histochem Cell Biol.* 2015; 143(6):565-74. (Biology)

## BD Biosciences

[bdbiosciences.com](http://bdbiosciences.com)

United States	Canada	Europe	Japan	Asia Pacific	Latin America/Caribbean
877.232.8995	866.979.9408	32.2.400.98.95	0120.8555.90	65.6861.0633	55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2015 BD

565798 Rev. 1

