

## Technical Data Sheet


**BB515 Rat Anti-Mouse CD195****Product Information**

<b>Material Number:</b>	<b>566208</b>
<b>Alternate Name:</b>	Ccr5; chemokine (C-C) receptor 5; C-C CKR-5; CC-CKR-5; CCR-5; AM4-7
<b>Size:</b>	25 µg
<b>Concentration:</b>	0.2 mg/ml
<b>Clone:</b>	C34-3448
<b>Immunogen:</b>	Mouse CCR5 aa. 9-30
<b>Isotype:</b>	Rat IgG2c, κ
<b>Reactivity:</b>	QC Testing: Mouse
<b>Storage Buffer:</b>	Aqueous buffered solution containing ≤0.09% sodium azide.

**Description**

The C34-3448 monoclonal antibody specifically binds to CD195 which is also known as, C-C chemokine receptor type 5 (CCR5). CD195 is a seven transmembrane-spanning G-protein-coupled receptor that belongs to the β-chemokine receptor family. CD195 regulates lymphocyte chemotaxis activation and transendothelial migration during inflammation. It signals in response to at least three chemokines: CCL3/MIP-1α, CCL4/MIP-1β, and CCL5/RANTES. CD195 is expressed on macrophages and some T-lymphocytes.

The antibody was conjugated to BD Horizon BB515 which is part of the BD Horizon Brilliant™ Blue family of dyes. With an Ex Max near 490 nm and an Em Max near 515 nm, BD Horizon BB515 can be excited by the blue laser (488 nm) laser and detected with a 530/30 nm filter. This dye has been exclusively developed by BD Biosciences and is up to seven times brighter than FITC with less spillover into the PE channel. Due to similar excitation and emission properties, BB515, FITC, and Alexa Fluor® 488 cannot be used simultaneously. It is not recommended to use BB515 in cocktails that include Streptavidin conjugates as it may cause high background.

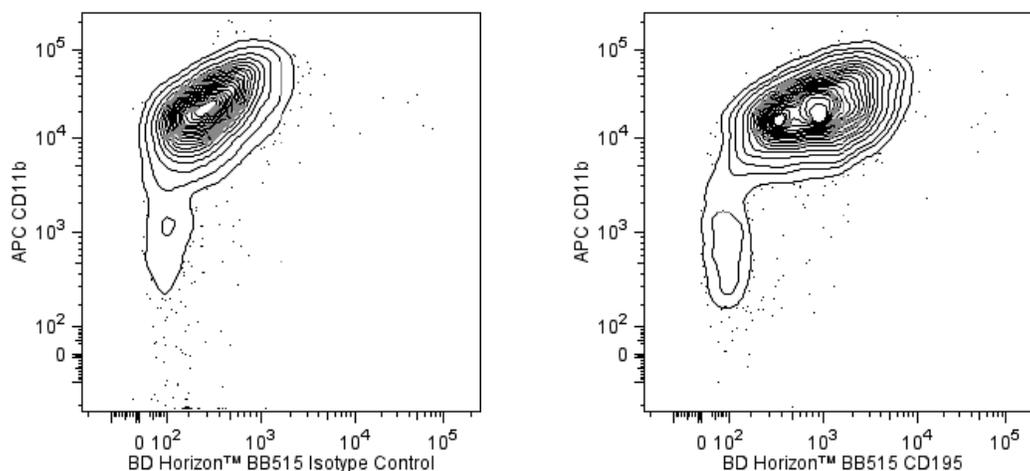

*Two-color flow cytometric analysis of CD195 expression on mouse peritoneal exudate cells. Mouse peritoneal exudate cells (PEC) were isolated 3 days post-stimulation by intraperitoneal injection of a 3% thioglycollate solution. The PEC were preincubated with Purified Rat Anti-Mouse CD16/CD32 antibody (Mouse BD Fc Block™) (Cat. No. 553141/553142). The cells were then stained with APC Rat Anti-Mouse CD11b (Cat. No. 553312/561690) and either BD Horizon™ BB515 Rat IgG2c, κ Isotype Control (Cat. No. 565091; Left Plot) or BD Horizon BB515 Rat Anti-Mouse CD195 (Cat. No. 565093/566208; Right Plot). Two-color flow cytometric contour plots showing the correlated expression of CD195 (or Ig Isotype control staining) versus CD11b were derived from gated events with the forward and side light-scatter characteristics of viable PEC. Flow cytometric analysis was performed using a BD™ LSR II Flow Cytometer System.*

**Preparation and Storage**

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ BB515 under optimum conditions and unconjugated antibody was removed.

**Application Notes****Application**

Flow cytometry	Routinely Tested
----------------	------------------

**BD Biosciences**

bdbiosciences.com

United States 877.232.8995 Canada 866.979.9408 Europe 32.2.400.98.95 Japan 0120.8555.90 Asia Pacific 65.6861.0633 Latin America/Caribbean 55.11.5185.9995

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

*Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.*

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.  
© 2017 BD. BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.



### Recommended Assay Procedure:

BD™ CompBeads can be used as surrogates to assess fluorescence spillover (Compensation). When fluorochrome conjugated antibodies are bound to CompBeads, they have spectral properties very similar to cells. However, for some fluorochromes there can be small differences in spectral emissions compared to cells, resulting in spillover values that differ when compared to biological controls. It is strongly recommended that when using a reagent for the first time, users compare the spillover on cells and CompBead to ensure that BD Comp beads are appropriate for your specific cellular application.

For optimal and reproducible results, BD Horizon Brilliant Stain Buffer should be used anytime two or more BD Horizon Brilliant dyes are used in the same experiment. Fluorescent dye interactions may cause staining artifacts which may affect data interpretation. The BD Horizon Brilliant Stain Buffer was designed to minimize these interactions. More information can be found in the Technical Data Sheet of the BD Horizon Brilliant Stain Buffer (Cat. No. 563794/566349) or the BD Horizon Brilliant Stain Buffer Plus (Cat. No. 566385).

For optimal results, it is recommended to perform 2 washes after staining with antibodies. Cells may be prepared, stained with antibodies and washed twice with wash buffer per established protocols for immunofluorescence staining, prior to acquisition on a flow cytometer. Performing fewer than the recommended wash steps may lead to increased spread of the negative population.

### Suggested Companion Products

Catalog Number	Name	Size	Clone
554656	Stain Buffer (FBS)	500 mL	(none)
554657	Stain Buffer (BSA)	500 mL	(none)
563794	Brilliant Stain Buffer	100 Tests	(none)
565093	BB515 Rat Anti-Mouse CD195	50 µg	C34-3448
565091	BB515 Rat IgG2c, κ Isotype Control	50 µg	A23-1
553141	Purified Rat Anti-Mouse CD16/CD32 (Mouse BD Fc Block™)	0.1 mg	2.4G2
553142	Purified Rat Anti-Mouse CD16/CD32 (Mouse BD Fc Block™)	0.5 mg	2.4G2
553312	APC Rat Anti-CD11b	0.1 mg	M1/70
561690	APC Rat Anti-CD11b	25 µg	M1/70
566349	Brilliant Stain Buffer	1000 Tests	(none)
566385	Brilliant Stain Buffer Plus	1000 Tests	(none)

### Product Notices

1. Since applications vary, each investigator should titrate the reagent to obtain optimal results.
2. An isotype control should be used at the same concentration as the antibody of interest.
3. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
4. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
5. The manufacture, use, sale, offer for sale, or import of this product is subject to one or more patents or pending applications. This product, and only in the amount purchased by buyer, may be used solely for buyer's own internal research, in a manner consistent with the accompanying product literature. No other right to use, sell or otherwise transfer (a) this product, or (b) its components is hereby granted expressly, by implication or by estoppel. Diagnostic uses require a separate license.
6. BD Horizon Brilliant Stain Buffer is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,575,303; 8,354,239.
7. Please refer to <http://regdocs.bd.com> to access safety data sheets (SDS).
8. Please refer to [www.bdbiosciences.com/us/s/resources](http://www.bdbiosciences.com/us/s/resources) for technical protocols.

### References

- Boring L, Gosling J, Monteclaro FS, Lulis AJ, Tsou CL, Charo IF. Molecular cloning and functional expression of murine JE (monocyte chemoattractant protein 1) and murine macrophage inflammatory protein 1α receptors: evidence for two closely linked C-C chemokine receptors on chromosome 9. *J Biol Chem*. 1996; 271(13):7551-7558. (Biology)
- Hayakawa Y, Smyth MJ. CD27 dissects mature NK cells into two subsets with distinct responsiveness and migratory capacity. *J Immunol*. 2006; 176(3):1517-1524. (Clone-specific: Flow cytometry)
- Kim SH, Cleary MM, Fox HS, Chantry D, Sarvetnick N. CCR4-bearing T cells participate in autoimmune diabetes. *J Clin Invest*. 2002; 110(11):1675-1686. (Clone-specific: Flow cytometry)
- McGuirk P, McCann C, Mills KH. Pathogen-specific T regulatory 1 cells induced in the respiratory tract by a bacterial molecule that stimulates interleukin 10 production by dendritic cells: a novel strategy for evasion of protective T helper type 1 responses by *Bordetella pertussis*. *J Exp Med*. 2002; 195(2):221-231. (Clone-specific: Flow cytometry)
- McLoughlin RM, Jenkins BJ, Grail D, et al. IL-6 trans-signaling via STAT3 directs T cell infiltration in acute inflammation. *Proc Natl Acad Sci U S A*. 2005; 102(27):9589-9594. (Clone-specific: Flow cytometry)
- Meyer A, Coyle AJ, Proudfoot AE, Wells TN, Power CA. Cloning and characterization of a novel murine macrophage inflammatory protein-1 alpha receptor. *J Biol Chem*. 1996; 271(24):14445-14451. (Biology)
- Nakae S, Iwakura Y, Suto H, Galli SJ. Phenotypic differences between Th1 and Th17 cells and negative regulation of Th1 cell differentiation by IL-17. *J Leukoc Biol*. 2007; 81(5):1258-1268. (Clone-specific: Flow cytometry)
- Napolitano M, Seamon KB, Leonard WJ. Identification of cell surface receptors for the Act-2 cytokine. *J Exp Med*. 1990; 172(1):285-289. (Biology)
- van Deventer HW, Wu QP, Bergstralh DT, et al. C-C chemokine receptor 5 on pulmonary fibrocytes facilitates migration and promotes metastasis via matrix metalloproteinase 9. *Am J Pathol*. 2008; 173(1):253-264. (Clone-specific: Flow cytometry)
- Zozulya AL, Reinke E, Baiu DC, Karman J, Sandor M, Fabry Z. Dendritic cell transmigration through brain microvessel endothelium is regulated by MIP-1α chemokine and matrix metalloproteinases. *J Immunol*. 2007; 178(1):520-529. (Clone-specific: Flow cytometry)