

Technical Data Sheet

BV510 Mouse Anti-Human CD19


Product Information

Material Number:	740164
Size:	50 µg
Clone:	HIB19
Alternative Name:	B4; B-lymphocyte antigen CD19; Leu-12
Reactivity:	Human (Tested in Development)
Isotype:	Mouse IgG1, κ
Application:	Flow cytometry (Qualified)
Concentration:	0.2 mg/ml
Workshop No.:	V CD19.11
Entrez Gene ID:	930
Storage Buffer:	Aqueous buffered solution containing ≤0.09% sodium azide.
Regulatory Status:	RUO

Description

The HIB19 monoclonal antibody specifically binds to the 95 kDa type I transmembrane CD19 glycoprotein. CD19 is expressed during all stages of B-cell maturation and differentiation, except on plasma cells. CD19 is also present on follicular dendritic cells. It is not found on T cells or on normal granulocytes. CD19 is a signal transduction molecule that regulates B cell development, activation, proliferation and differentiation. It associates with the complement receptor 2 (CD21), TAPA-1 (CD81), Leu 13, and/or MHC class II to form a signal transduction complex on the surface of B cells. Anti-CD19 clone HIB19 partially blocks the binding of clone B43, another CD19-specific monoclonal antibody.

The antibody was conjugated to BD Horizon™ BV510 which is part of the BD Horizon Brilliant™ Violet family of dyes. With an Ex Max of 405-nm and Em Max at 510-nm, BD Horizon BV510 can be excited by the violet laser and detected in the BD Horizon V500 (525/50-nm) filter set. BD Horizon BV510 conjugates are useful for the detection of dim markers off the violet laser.

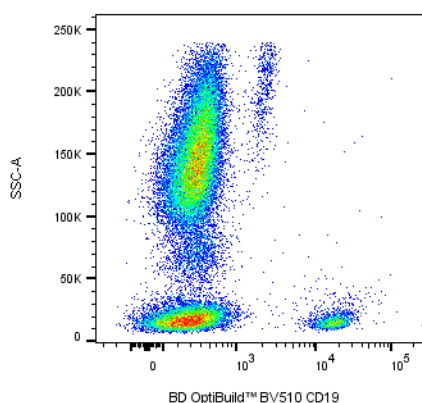

Multiparameter flow cytometric analysis using BD OptiBuild™ BV510 Mouse Anti-Human CD19 antibody (Cat. No. 740164) on human peripheral lymphocytes. Flow cytometry was performed using a BD FACSCelesta™ Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze. The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography. The antibody was conjugated with BD Horizon BV510 under optimal conditions that minimize unconjugated dye and antibody.

Recommended Assay Procedure

For optimal and reproducible results, BD Horizon Brilliant Stain Buffer should be used anytime two or more BD Horizon Brilliant dyes (including BD OptiBuild Brilliant reagents) are used in the same experiment. Fluorescent dye interactions may cause staining artifacts which may affect data interpretation. The BD Horizon Brilliant Stain Buffer was designed to

minimize these interactions. More information can be found in the Technical Data Sheet of the BD Horizon Brilliant Stain Buffer (Cat. No. 563794).

Suggested Companion Products

Catalog Number	Name	Size	Clone
562946	BV510 Mouse IgG1, k Isotype Control	50 µg	X40
554656	Stain Buffer (FBS)	500 mL	
554657	Stain Buffer (BSA)	500 mL	
563794	Brilliant Stain Buffer	100 Tests	
555899	Lysing Buffer	100 mL	
349202	Lysing Solution 10X Concentrate	100 NA	
564219	Human BD Fc Block™	50 mg	

Product Notices

1. This antibody was developed for use in flow cytometry.
2. The production process underwent stringent testing and validation to assure that it generates a high-quality conjugate with consistent performance and specific binding activity. However, verification testing has not been performed on all conjugate lots.
3. Researchers should determine the optimal concentration of this reagent for their individual applications.
4. An isotype control should be used at the same concentration as the antibody of interest.
5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
7. Please refer to www.bdbiosciences.com/us/s/resources for technical protocols.
8. BD Horizon Brilliant Stain Buffer is covered by one or more of the following US patents: 8,110,673; 8,158,444; 8,575,303; 8,354,239.
9. BD Horizon Brilliant Violet 510 is covered by one or more of the following US patents: 8,575,303; 8,354,239.

References

Bradbury LE, Goldmacher VS, Tedder TF. The CD19 signal transduction complex of B lymphocytes. Deletion of the CD19 cytoplasmic domain alters signal transduction but not complex formation with TAPA-1 and Leu 13. *J Immunol.* 1993; 151(6):2915-2927. (Biology: Flow cytometry).

Favaloro EJ, Moraitis N, Koutts J, Exner T, Bradstock KF. Endothelial cells and normal circulating haemopoietic cells share a number of surface antigens. *Thromb Haemost.* 1989; 61(2):217-224. (Biology: Flow cytometry).

Knapp W. W. Knapp .. et al., ed. *Leucocyte typing IV : white cell differentiation antigens.* Oxford New York: Oxford University Press; 1989; :1-1182.

Nadler LM, Anderson KC, Marti G, et al. B4, a human B lymphocyte-associated antigen expressed on normal, mitogen-activated, and malignant B lymphocytes. *J Immunol.* 1983; 131(1):244-250. (Biology: Flow cytometry).

Schlossman SF, Stuart F, Schlossman .. et al., ed. *Leucocyte typing V : white cell differentiation antigens : proceedings of the fifth international workshop and conference held in Boston, USA, 3-7 November, 1993.* Oxford: Oxford University Press; 1995; .

Uckun FM, Muraguchi A, Ledbetter JA, et al. Biphenotypic leukemic lymphocyte precursors in CD2+CD19+ acute lymphoblastic leukemia and their putative normal counterparts in human fetal hematopoietic tissues. *Blood.* 1989; 73(4):1000-1015. (Biology: Flow cytometry).

BD Biosciences

bdbiosciences.com

United States
877.232.8995

Canada
888.268.5430

Europe
32.53.720.550

Japan
0120.8555.90

Asia Pacific
65.6861.0633

Latin America/Caribbean
0800.771.7157



For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for a patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton Dickinson and Company is strictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

©2020 BD. All rights reserved. Unless otherwise noted, BD, the BD Logo and all other trademarks are the property of Becton, Dickinson and Company or its affiliates.